

photoLab S12



Operating Instructions

Part 1: General Information

Part 2: Functional Description

Contents

1	Photometers	2
1.1	Photometry	2
1.2	The Photometers	3
2	Photometric Test Kits	3
2.1	Basic Principle	3
2.1.1	Spectroquant® Cell Tests	4
2.1.2	Spectroquant® Reagent Tests	4
2.2	Notes for Practical Use	5
2.2.1	Measuring Range	5
2.2.2	Influence of pH	7
2.2.3	Influence of Temperature	7
2.2.4	Time Stability	7
2.2.5	Influence of Foreign Substances	8
2.2.6	Dosing of Reagents	8
2.2.7	Shelf-life of the Reagents	9
3	Sample Preparation	9
3.1	Taking Samples	9
3.2	Preliminary Tests	10
3.3	Dilution	10
3.4	Filtration	11
3.5	Homogenization	12
3.6	Decomposition	12
4	Pipetting System	14
5	Analytical Quality Assurance (AQA)	15
5.1	Quality Control at the Manufacturer	15
5.2	Quality Control for the User	16
5.2.1	Checking the Photometer	17
5.2.2	Checking the Overall System	17
5.2.3	Checking the Pipettes	18
5.2.4	Checking Thermoreactors	18
5.2.5	Testing for Handling Errors	19
5.3	Determination of Sample Influences	19
5.4	Definition of Errors	20

1 Photometers

1.1 Photometry

When a beam of light is transmitted through a coloured solution, then this beam loses its intensity, in other words a part of the light is absorbed by the solution. Depending on the substance in question, this absorption occurs at specific wavelengths.


Monochromators (e.g. narrow-band interference filters, lattices) are used to select the wavelength from the total spectrum of a tungsten-halogen lamp (VIS spectrum), a deuterium lamp (UV spectrum) or, respectively, a xenon lamp.

The intensity of the absorption can be characterized using the transmittance T (or, respectively, T in percent).

$$T = I/I_0$$

I_0 = Initial intensity of the light

I = Intensity of the transmitted light

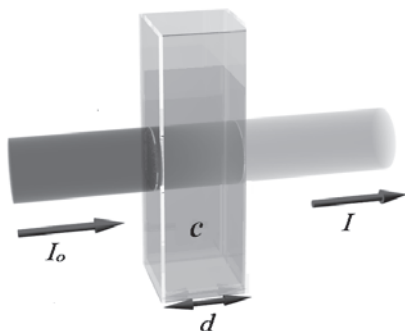

If the light is not absorbed at all by a solution, then this solution has a transmittance of 100 %; a complete absorption of the light in the solution means 0 % transmittance.

The measure generally used for the absorption of light is the absorbance (A), since this correlates directly with the concentration of the absorbing substance. The following connection exists between absorbance and transmittance:

$$A = -\log T$$

Experiments by BOUGUER (1698–1758) and LAMBERT (1728–1777) showed that the absorbance is dependent on the thickness of the absorbing layer of the cell used. The relationship between the absorbance and the concentration of the analyte in question was discovered by BEER (1825–1863). The combination of these two natural laws led to the derivation of *Lambert-Beer's law*, which can be described in the form of the following equation:

$$A = \epsilon_{\lambda} \cdot c \cdot d$$

ϵ_{λ} = Molar absorptivity, in $l/mol \times cm$

d = Path length of the cell, in cm

c = Concentration of the analyte, in mol/l

1 Photometers

1.2 The Photometers

The photometers that belong to the Spectroquant® Analysis System differ from conventional photometers in the following important aspects:

- The calibration functions of all test kits are electronically stored.
- The measurement value can be immediately read off from the display in the desired form.
- The method for the test kits (Cell Tests **and** reagent tests) belonging to the Spectroquant® analysis system is automatically selected via the scanning of the bar code.
- All cells formats used are automatically identified and the correct measuring range is selected automatically.
- Instrument-supported AQA ensures that measurement results can be used as secure, reproducible, and recognized analytical results.
- New methods can be downloaded from the internet site <http://photometry.merck.de> and permanently stored in the instrument.

For technical data and instructions for use please refer to the section “Function description” or can also be found on the internet.

2 Photometric Test Kits


2.1 Basic Principle

By means of reagents, the component of a sample to be analyzed is converted into a coloured compound in a specific reaction. The reagents or reagent mixtures contain – in addition to the reagent selective for a parameter to be determined – a number of auxiliary substances that are essential for the course of the reaction. These include, for example, buffers for adjusting the pH to the optimal value for the reaction, and masking agents that suppress or minimize the influence of interfering ions.

The colour reactions are in most cases based on standardized analytical methods specifically optimized in terms of ease of use, a low working effort, and shorter reaction times. Furthermore, methods cited in the literature or developed by ourselves are also used. Details on the respective reference procedures are stated in the package insert or else in the parameter overview.

2 Photometric Test Kits

2.1.1 Spectroquant® Cell Tests


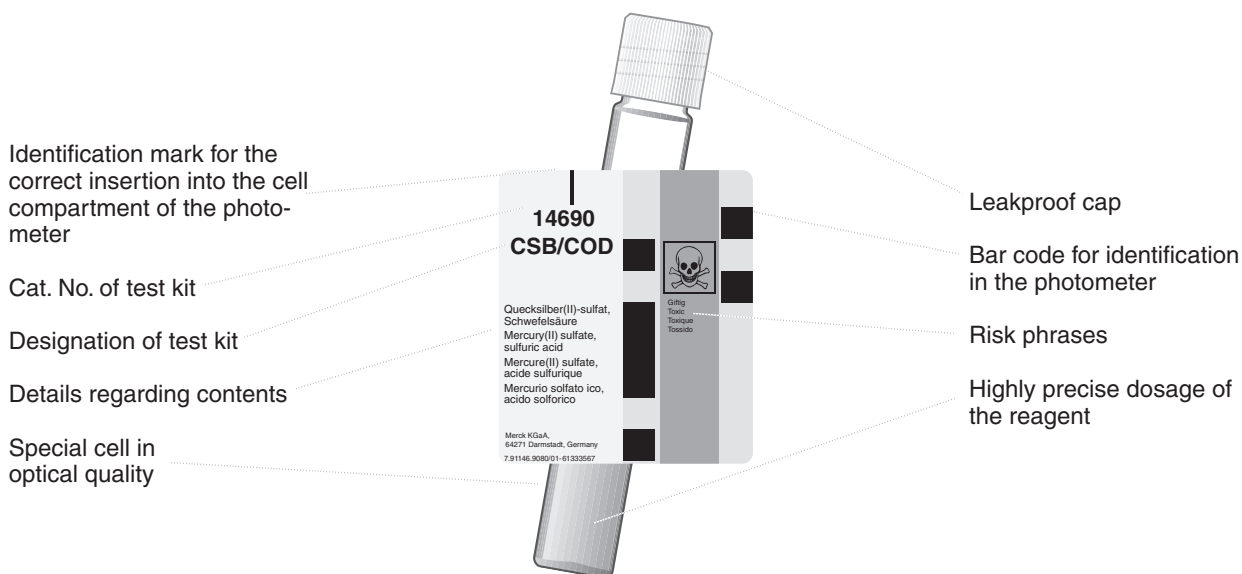


Additional reagent(s)

Certain cell tests, e. g. COD or nitrite, already contain all necessary reagents in the cells, and the sample must merely be added with a pipette. In other tests, however for reasons of chemical compatibility it is necessary to separate the test into two or three different reagent mixtures. In such cases, besides the sample a metered reagent must also be added.

2.1.2 Spectroquant® Reagent Tests

The principle behind the reagent tests is that the reagents necessary for the colour reaction are combined in the form of liquid concentrates or solid-substance mixtures. A few drops of the reagent concentrate are added to the sample. This means that there is no need to dilute the sample, which in turn enhances the sensitivity of the detection. The procedure generally used in classical photometry by which the sample is made up to a defined volume in a volumetric flask is dispensed with.

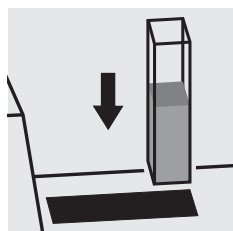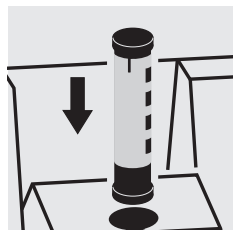

The method is selected automatically by means of the scanning of the bar code by the AutoSelector.

All cells formats used are automatically identified and the correct measuring range is selected automatically.

Subsequently the result is automatically shown on the display.

2 Photometric Test Kits

2.2 Notes for Practicle Use


2.2.1 Measuring range

The intensity of the colour of a solution, measured as the absorbance, is proportional to the concentration of the respective analyte only within a specific range. This measuring range (effective range) is electronically stored in the photometers for each individual test kit .

Below the specified measuring range, either a different cell or else another procedure must be used. The **lower limit of the measuring range** either takes the form of nonlinearity of the calibration curve, as shown in the figure, or else is given by the method detection limit. The **method detection limit** of an analytical method is the lowest concentration of the analyte in question that can be measured quantitatively with a defined degree of probability (e.g. 99 %).

The **upper limit of the measuring range** is the point at which the linear correlation between the concentration and the absorbance ends. In such a case the sample must be diluted accordingly so that it lies ideally in the middle of the effective range (least-error measurement).

In photometry it is conventional practice to measure against the reagent blank value. Here the analysis is carried out "blind", i.e. without any analyte added. Instead of the sample volume, the corresponding quantity of distilled or DI water is used. This **reagent blank value is prestored** in the photometers belonging to the Spectroquant® analysis system, which means that - due to the high batch reproducibility - it is possible to dispense with a separate measurement of the reagent blank. At the lower limit of the measuring range, the accuracy of the determination can be enhanced by performing the measurement against a separately prepared reagent blank.

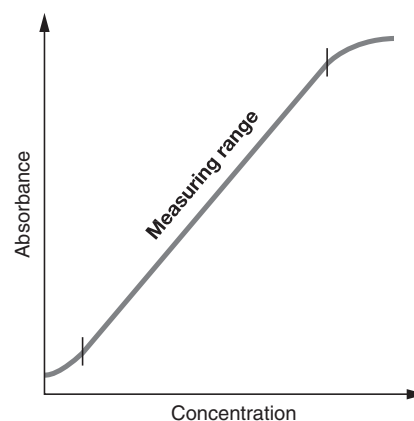

2 Photometric Test Kits

In some cases the intensity of the colour of the solution and thus the absorbance can drop again when **very high concentrations of the analyte** are present. The examples are listed in the following table. The values indicated in the display are correct up to the concentrations specified in the third column, and false measuring values are obtained above these concentrations. In such a case it is necessary to conduct a plausibility check by running preliminary tests using test strips or dilution.

Cat. No.	Method	Correct indication of result up to sample conc.	Farbänderung
14739	Ammonium CT	100 mg/l	turquoise instead of green
14558	Ammonium CT	500 mg/l	turquoise instead of green
14544	Ammonium CT	1000 mg/l	turquoise instead of green
14559	Ammonium CT	5000 mg/l	turquoise instead of green
14752	Ammonium Test	100 mg/l	turquoise instead of green
00683	Ammonium Test	2500 mg/l	turquoise instead of green
00605	Bromine Test	50 mg/l	yellow instead of red
00595	Chlorine CT	25 mg/l	yellow instead of red
00597	Chlorine CT	25 mg/l	yellow instead of red
00598	Chlorine Test	25 mg/l	yellow instead of red
00602	Chlorine Test	25 mg/l	yellow instead of red
00599	Chlorine Test	25 mg/l	yellow instead of red
00086/ 87/88	Chlorine Test	300 mg/l	yellow instead of red
00608	Chlorine Dioxide Test	15 mg/l	yellow instead of red
14553	Copper CT	50 mg/l	turquoise instead of blue
14767	Copper Test	50 mg/l	turquoise instead of blue
14557	Fluoride CT	5 mg/l	brown instead of violet
14598	Fluoride Test	5 (50) mg/l	brown instead of violet
00606	Iodine Test	50 mg/l	yellow instead of red
01632	Monochloramine Test	300 mg/l	turquoise instead of green
00607	Ozone Test	15 mg/l	yellow instead of red
14551	Phenol CT	100 mg/l	weakening of colour
14831	Silver Test	5 mg/l	no change (flocculation)

2 Photometric Test Kits

2.2.2 Influence of pH

Chemical reactions follow an optimal course only within a certain pH range. The reagents contained in the test kits produce an adequate buffering of the sample solutions and ensure that the pH optimal for the reaction in question is obtained.

Strongly acidic ($\text{pH} < 2$) and strongly alkaline ($\text{pH} > 12$) sample solutions can prevent the pH from being adjusted to an optimal range, since under certain circumstances the buffering capacity of the test-kit reagents may not be sufficient. Any necessary correction is made by the dropwise addition of diluted acid (reduces the pH) or diluted lye (raises the pH), testing the pH with suitable indicator strips after each drop is added. The addition of the acid or lye results in a dilution of the test solution. When up to five drops are added to 10 ml of sample, the change in the volume can be neglected, since the resultant error is lower than 2 %. The addition of larger quantities should be duly considered by adjusting the sample volume accordingly.


The specified pH values for the sample solution and, wherever applicable, for the measurement solution are defined in the respective package inserts and in the analysis instructions in chapter 3 of the manual.

2.2.3 Influence of Temperature

The temperature of the sample solution and the reagents may have an effect on the colour reaction and thus on the measurement result. The typical temperature course is illustrated in the figure.

If the sample temperature is lower than 15 °C, false-low results must be reckoned with. Temperatures exceeding 30 °C generally influence the stability of the compound that is formed in the reaction. The optimal temperature for the colour reaction is stated in the package inserts of the respective Spectroquant® test kits.


Attention! After thermic decomposition procedures, the determination of COD or total contents of nitrogen, phosphorus, or metal, a sufficient waiting time must be allowed for to permit the solution cool to room temperature.

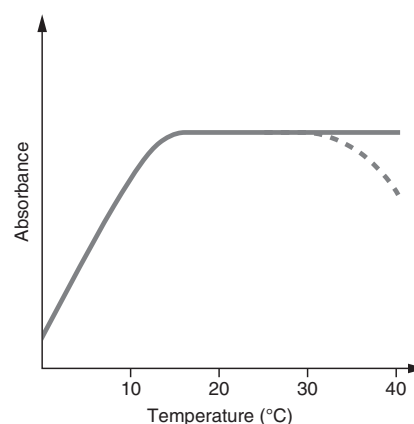

2.2.4 Time Stability

Most of the colour reactions require a certain time to reach the maximum colour intensity. The solid curve in the figure at the right gives a schematic impression of a typical time course. The behaviour of relatively instable colour reactions with time is shown by the dotted curve.

The reaction time specified in the working instructions refers to the period of time from the addition of the last reagent until the actual measurement. In addition, the package inserts for the individual test kits also state the time interval in which the measurement value does not change. The maximum time interval is 60 minutes; this time should not be exceeded, even in the case of stable colour reactions.

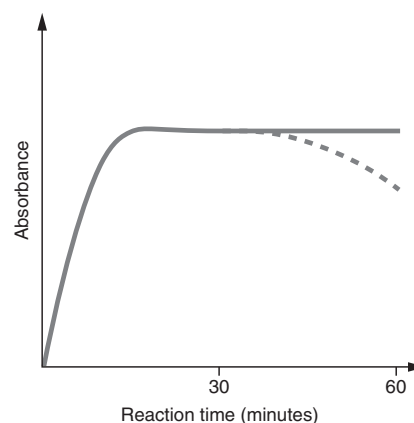

2 Photometric Test Kits

2.2.5 Influence of Foreign Substances

Foreign substances in the sample solution can

- raise the measurement value as a result of an amplification of the reaction
- lower the measurement value as a result of a prevention of the reaction.

A quantification of this effects is stated in tabular form in the respective package inserts for the most important foreign ions. The tolerance limits have been determined for the individual ions; they may not be evaluated cumulatively.

Suitability for use in seawater

A tabular survey (see appendix 1) provides information on the suitability of the tests in connection with seawater and also on the tolerances for salt concentrations.


2.2.6 Dosing the Reagents



Small amounts of liquids are dosed by counting the number of drops from a leakproof bottle.



When using dropper bottles it is extremely important that the bottle be held vertically and that the drops be added slowly (approx. 1 drop per second). If this is not observed, the correct drop size and thus the correct amount of reagent are not achieved.

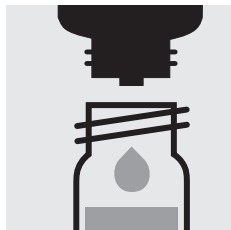


A positive-displacement pipette should be used for larger quantities of liquid or for the exact dosage of smaller reagent quantities. In these cases the reagent bottles are not fitted with a dropper insert.



Solid substances are dosed either with the dose-metering cap or with microspoons that are integrated into the screw cap of the respective reagent bottle. The dose-metering cap is used for solid reagents or reagent mixtures that are free-flowing.

In all other cases the substances are dosed with the microspoon.

In this case it is necessary to add only level microspoonfuls. To this end the spoon must be drawn over the brim of the reagent bottle.

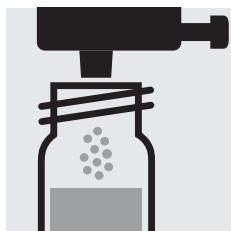

At the first use replace the black screw cap of the reagent bottle by the dose-metering cap.

Hold the reagent bottle vertically and, at each dosage, press the slide all the way into the dose-metering cap. Before each dosage ensure that the slide is completely retracted.



Reclose the reagent bottle with the black screw cap at the end of the measurement series, since the function of the reagent is impaired by the absorption of atmospheric moisture.

2 Photometric TestKits

2.2.7 Shelf-life of the Reagents

The Spectroquant® test kits are in most cases stable for 3 years when stored in a cool, dry place. A few test kits have a lower shelf-life of 18 or 24 months or must else be stored in a refrigerator.

COD Cell Tests must be stored protected from light.

The expiry date of the package unit is printed on the outer label. The shelf-life may become reduced when the reagent bottles are not reclosed tightly after use or when the test kit is stored at temperatures higher than those specified.

3 Sample Preparation

Sample preparation covers all the steps necessary before the actual analysis can be performed.

3.1 Taking Samples

The taking of samples is the first and **most important step** on the way to obtaining the correct analysis result. Not even the most exact method of analysis can correct any mistakes made in the taking of the sample. The objective of the sampling procedure is to gain a sample with a representative composition. The most important precondition for **gaining a representative sample** is the identification of the suitable sampling site. Here it must be borne in mind that the solution to be investigated can display varying concentrations in different places at different times.

In sampling, a distinction is made between manual and automatic methods. In many cases a true picture of the average composition of the sample can be obtained only once several individual samples have been collected; this can be done manually or with an automatic sampler.

Clean plastic containers with a volume of 500 or 1000 ml are suitable for collecting samples. They should be rinsed several times, under vigorously shaken, with the water to be investigated, and then filled free of air bubbles and immediately closed tightly. The containers must be protected against the effects of air and heat and then be forwarded for the further analytical steps as soon as possible. In exceptional cases, preservation measures in the form of short-term refrigeration at +2 to +5 °C and chemical conservation can be taken.

Parameter	Preservation
COD	+2 to +5 °C max. 24 h or -18 °C max. 14 days
N compounds: NH ₄ -N, NO ₃ -N, NO ₂ -N	analyze immediately, only in exceptional case +2 to +5 °C max. 6 h
P compounds: PO ₄ -P, P total	short-term storage, no preservation; with nitric acid to pH 1, max. 4 weeks
Heavy metals	short-term storage, no preservation; with nitric acid to pH 1, max. 4 weeks

3 Sample Preparation

3.2 Preliminary Tests

Correct measurement results can be obtained only within the measuring range specified for each individual parameter. When dealing with sample solutions of an unknown concentration, it is advisable to establish whether the sample concentration is indeed within the specified measuring range, ideally roughly in the middle of the range.

Preliminary tests enhance the analytical reliability and make the determination of the necessary dilution ratios in the case of high concentrations easier. **Merckoquant® Test Strips** are very well suited for preliminary tests.

3.3 Dilution

Dilution of samples is necessary for two reasons:

- The concentration of the parameter under investigation is too high, i. e. it lies outside the measuring range.
- Other substances contained in the sample interfere with the determination (matrix interference); false-high or false-low results may ensue.

The following auxiliaries are absolute prerequisites for the dilution of the sample:

- Volumetric flasks of varying sizes (e. g. 50, 100 and 200 ml)
- Positive-displacement pipette
- Distilled or DI water.

Only dilutions carried out with these auxiliary products are of sufficient reliability in the area of trace analysis, to which photometry belongs (for the simplified procedure see page 11).

An important aspect here is that once the volumetric flask has been filled up to the mark with distilled water the flask is closed and the contents are thoroughly mixed.

The **dilution factor (D_F)** resulting from the dilution procedure is calculated as follows:

$$D_F = \frac{\text{Final volume (total volume)}}{\text{Initial volume (sample volume)}}$$

The analytical result is subsequently multiplied by the dilution factor.

A calculation can be dispensed with when the dilution is programmed into the photometer. The **dilution number** (see the table on page 11) is entered and the measurement value is subsequently calculated correctly and immediately displayed.

3 Sample Preparation

All dilutions should be made in such a way that the measurement value lies in the middle of the measuring range. As a rule, the dilution factor should never be higher than 100. In the event that yet higher dilutions become necessary all the same, then this must be done in two separate steps.

Example

Step 1: Make up 2 ml of sample to 200 ml with distilled water;
 $D_F = 100$, dilution number 1+99

Step 2: Take 5 ml of the above solution and make up to 100 ml;
 $D_F = 20$, dilution number 1+19

The dilution factor for the total dilution is calculated by multiplying the individual dilutions:

$$D_{F_{\text{total}}} = D_{F_1} \times D_{F_2} = 100 \times 20 = 2000, \text{ dilution number } 1+1999$$

Simplified procedure

Dilutions up to 1:10 can also be prepared without volumetric flasks in a glass beaker, measuring the volumes of the sample and the dilution water using a previously calibrated positive-displacement pipette (see table for instructions).

Desired dilution	Volume of sample [ml]	Volume of dist. water [ml]	Dilution factor	Dilution number
1:2	5	5	2	1+1
1:3	5	10	3	1+2
1:4	2	6	4	1+3
1:5	2	8	5	1+4
1:10	1	9	10	1+9

3.4 Filtration

Strongly turbid samples require pretreatment before they can be determined in a photometer, since the effect of turbidity can result in considerable variations in the measurement values and in false-high readings. Care must be taken here to ensure that the substance to be determined is not contained in the suspended material, in which case a sample decomposition must be carried out.


Compounds that always occur in dissolved form (for example ammonium, nitrate, nitrite, chlorine, chloride, cyanide, fluoride, orthophosphate, and sulfate) permit a previous filtration, even when the sample solution is strongly turbid.

Weak turbidity is eliminated by the **automatic turbidity-correction** feature built into the photometer (see Function description, "Device set-up/Correction function"); in such cases it is not necessary to filter the sample before analysis.


As a measure to distinguish between dissolved and undissolved water-borne substances, the water sample can be filtered through a simple paper filter. Following the recommendations stated in the reference methods, membrane filters with a pore size of 0.45 µm are required for fine filtration.

3 Sample Preparation


Procedure for microfiltration

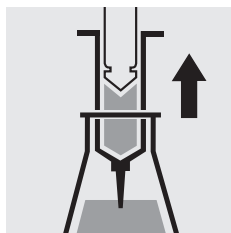


Draw out the liquid to be filtered with the syringe.

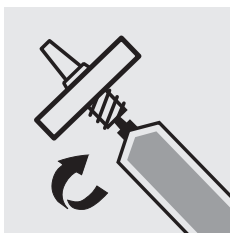

Screw the syringe tightly into the front side of the membrane-filter attachment.

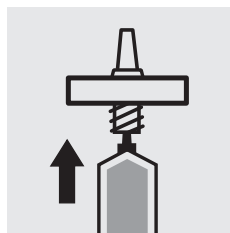

Hold the syringe upright and slowly depress the piston upwards until the membrane-filter is fully wetted free of air bubbles.

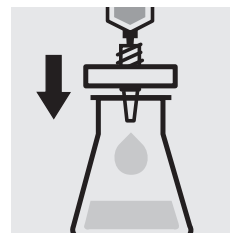


Filter the contents of the syringe into the intended glass vessel.

3.5 Homogenization

As a measure to ensure that a representative sample can be taken in the presence of suspended matter in the water sample in question, for certain parameters - e.g. COD and the total content of heavy metals - the sample must be homogenized. This must be carried out using a high-speed blender (2 minutes at 5000 – 20 000 rpm and taking the sample while stirring).

3.6 Decomposition

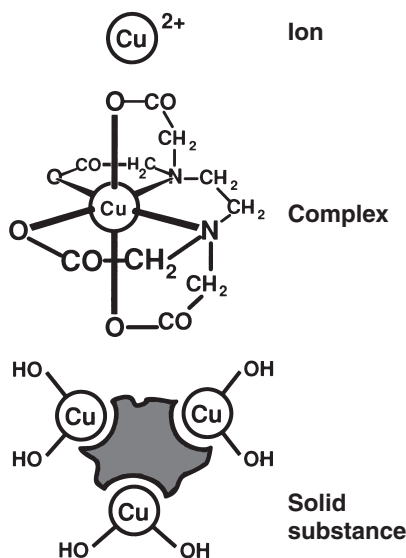


Water-borne substances can be present in the sample for investigation in a variety of forms: as the ion, bound more or less solidly in a complex, or as a solid substance.

3 Sample Preparation

The manner in which the sample is pretreated enables the three proportions to be distinguished from each other. This can be illustrated using a copper-containing wastewater sample as an example.

Example

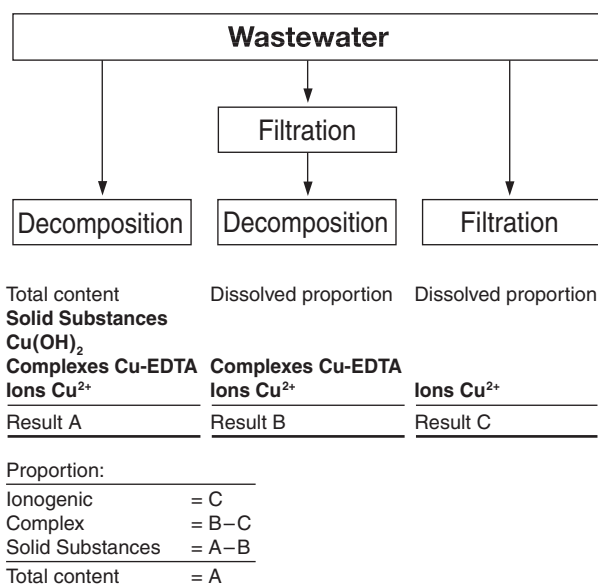

Decomposition converts the substance to be determined into an analyzable form. In most cases, decomposition agents take the form of acids in combination with oxidizing agents; in exceptional cases (e. g. in the determination of total nitrogen) an alkaline decomposition is more effective. The type of decomposition procedure used depends on the analyte to be determined and the sample matrix.

The ready-to-use sample-decomposition products **Spectroquant® Crack Set 10** and **20** are suited for the preparation of the sample materials for the determinations stated in the table.

The decomposition processes are carried out in the **Spectroquant® thermoreactor** (capacity: 12 or 24 decomposition cells) at 120 °C or, respectively, 100 °C. Details regarding the heating times and further treatment can be found in the package inserts contained in the **Spectroquant® Crack Set** packs.

Determination of	Sample preparation with
Total phosphorus*	Crack Set 10 / 10C**
Total chromium* [= sum of chromate and chromium(III)]	Crack Set 10 / 10C
Total metal [= sum of free and complex-bound metal]	Crack Set 10 / 10C
Total nitrogen*	Crack Set 20


* The decomposition reagents are already contained in the packs of the respective cell tests.

** Decomposition cells are included in the pack; empty cells are required for the decomposition for Crack Sets 10 and 20.

3 Sample Preparation

In the event that the sample to be analyzed is a highly contaminated material (high proportion of organic substances) or water-insoluble samples, decomposition using concentrated acids and other agents is indispensable. Corresponding examples from the **collection of applications** for real samples are available on request.

The necessity for decomposition can be checked according to the following diagram:

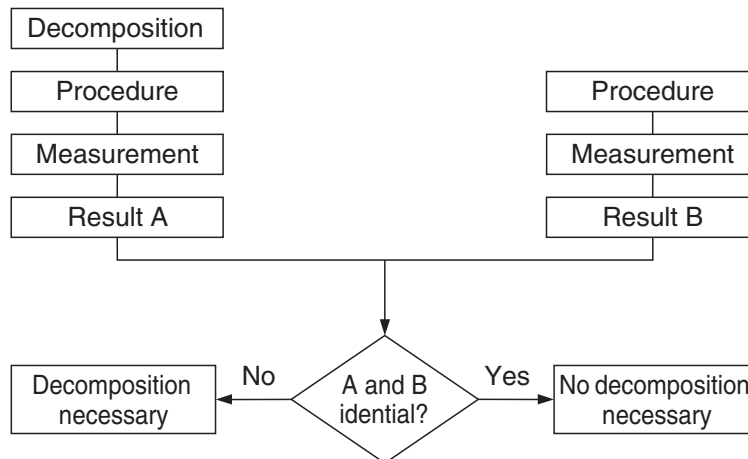

For wastewater with a consistent composition, this check as a rule need be carried out only once. It is, however, advisable to check the result periodically.

4 Pipetting System

Positive-displacement pipettes permit

- an exact dosage of the sample volume
- a precise measurement of sample and reagent volumes and of the volumes of water for dilution purposes.

Pipettes of varying volumes and also ones with a fixed volume are available.

Sources of error and hints on how to avoid them:

- Closely follow the instructions for use contained with the pipette in question.
- Check the pipetted volumes
 - a) by weighing using analytical scales (weighing accuracy ± 1 mg), 1 ml of water at 20 °C = 1.000 g ± 1 mg
 - b) using Spectroquant® PipeCheck; this is a photometric check of the pipette, and scales are not necessary (see section "AQA").
- Avoidance of spread effects by rinsing the pipette several times with the solution to be pipetted.
- Always exchange the pipette tip.
- Draw up the liquid slowly and depress piston completely to discharge the liquid.

5 Analytical Quality Assurance (AQA)

The objective of analysis must always be to determine the true content of the analyte in question as accurately and precisely as possible.

Analytical Quality Assurance represents a suitable and indispensable method by which the quality of the user's own work can be assessed, errors in the measurement system diagnosed, and the comparability with the results obtained using the respective reference methods demonstrated.

Details regarding the necessity of AQA can be found in the in Memorandum A 704 of the German Association for the Water Sector, Wastewater, and Waste Materials (Deutsche Vereinigung für Wasserwirtschaft, Abwasser und Abfall e.V., DWA) and in the corresponding self-control/self-monitoring regulations of the German federal states (available in english).

Causes for errors can include:

- the working materials used
- the handling
- the sample under investigation.

These errors have effects on both the accuracy and precision of the results obtained.


5.1 Quality Control at the Manufacturer

Photometers and photometric test kits possess specifications that are adhered to and above all else also documented by the manufacturer.

The **certificate for the photometer** enclosed with each device documents the quality of the measuring device.

Prüfprotokoll Test record		Seite 1 von 1 page 1 of 1
Photometer / photometer	Modell / Model Serien-Nr. / Serial no.	SQ NOVA 60 Akku 07046913
Transponderfunktion / funktion of transponder		
Korrektles Erlesen eines Testdatensatzes / correct reading of test data set		ok
Selbsttest / Self Check		
Signalstärk ohne Küvette / signal adjustment without cell		ok
LS-Check / LS-Check		
Korrekte Erkennung von Test-Barcodes / correct identification of test barcodes		ok
Nullabgleich Rechteckküvette 10 mm / Zero adjustment rectangular cell 10 mm		ok
Nullabgleich Rechteckküvette 20 mm / Zero adjustment rectangular cell 20 mm		ok
Nullabgleich Rechteckküvette 50 mm / Zero adjustment rectangular cell 50 mm		ok
Nullabgleich Rundküvette / Zero adjustment round cell		ok
Bei allen Filterpositionen Abgleich auf Extinktion E = 0 mit unchloriertem Wasser als Messlösung / at all filterpositions adjustment to absorbance A = 0 using dechlorinated water as measurement solution		ok
Photometrische Richtigkeit / photometric accuracy		
Extinktion E einer Testlösung in Rundküvette / absorbance A of test solution in round cell		
Wellenlänge / wavelength (nm): 695	Softwert E / nominal value A: 0.385	
Toleranz E / toleranz A: ±0.020	Messwert E / measured value A: 0.379	ok
Linearität / linearity		
Extinktionsdaten von 2-Platzfiltern in separater (E1-E2) und kombinierter Anordnung (E12) / absorbance data of 2-place filters in separate (E1,E2) and combined configuration (E12)		
Messwerte E / Measured values A (AASnm):	E1 = 0.998 E2 = 1.015 E12 = 2.001	
Anforderung / requirement:	0.02 ≤ (E1 + E2 - E12) ≤ 0.02	ok
Elektrische Sicherheit nach IEC 1010 / electrical safety according to IEC 1010		
Keine visuellen Mängel, keine Girax, keine losen Teile und Befestigungen / no visual flaws, no burrs, no loose parts and fastenings		ok
Datum / Date:	30.01.2007	Prüfer / Tester: Michael Dobry
Dieses Dokument wurde mittels EDV erstellt und ist ohne Unterschrift gültig / This document has been generated using electronic data processing and is valid without signature		

5 Analytical Quality Assurance (AQA)

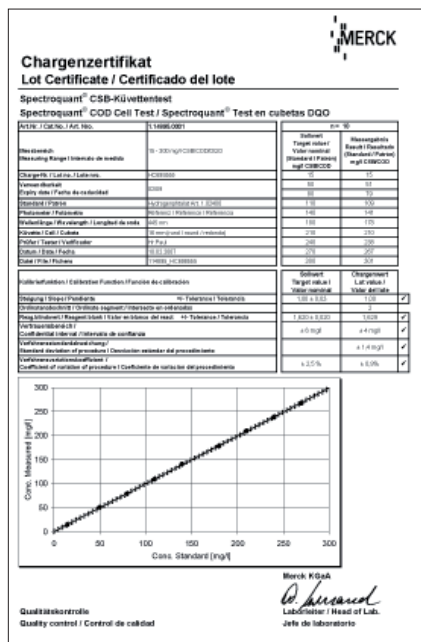

The **certificate for the test kit**, available for each lot produced, documents the quality of the reagents contained in the test kit.

Calibration function:
The calculated function must agree, within specified tolerances, with the function electronically stored in the photometer.

Confidence interval:
Maximum deviation from the desired value over the entire measuring range; every measurement value can be affected by this deviation; this parameter is a measure for the accuracy.

Standard deviation for the procedure:
Measurement for the dispersion of the measurement values over the entire measuring range, expressed in ± mg/l.


Coefficient of variation for the procedure:
Measurement for the dispersion of the measurement values over the entire measuring range, expressed in %. The smaller the standard deviation/coefficient of variation for the procedure, the more pronounced the linearity of the calibration curve.

5.2 Quality Control for the User

A complete check comprises the entire system, i. e. the working equipment and the mode of operation. The photometer offers an optimum degree of support in this regard, in the form of the different quality mode. The instrument, or the whole system (including reagents and all accessories) will be checked, depending on which quality mode selected. All of checking operations can thus be supported by the photometer and the check values accordingly documented as per GLP (Good Laboratory Practice) recommendations (see Function description, "Analytical Quality Assurance").


The following diagram provides an overview regarding internal quality-assurance aspects:

Checking the working equipment

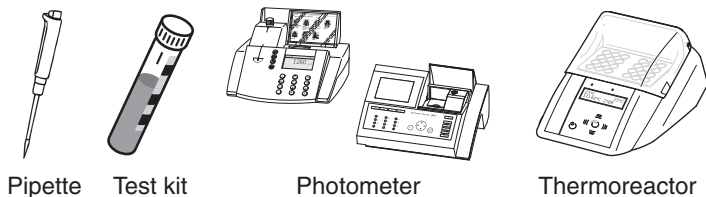

= Test for the overall system

Checking the handling operations

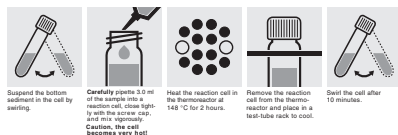

Influence of the sample

Test for recovery

5 Analytical Quality Assurance (AQA)

5.2.1 Checking the Photometer


As soon as the photometer is activated it is running a Self-Check. This means the hardware and the software of the photometer is checked and compared with internal standards.

As soon as the photometer is activated it is running a Self-Check. This means the hardware and the software of the photometer is checked and compared with internal standards.

The photometer itself is checked in the **AQA 1 mode** with the **Spectroquant® PhotoCheck**: the pack includes round cells containing stable test solutions (**secondary standards**) for checking the photometer at the 445, 525, and 690 nm wavelengths. The test solutions are measured in a **reference photometer** monitored with **primary standards**, and the certificate stating the absorbance values is enclosed with the package unit. These desired values with the permissible tolerances are entered into the photometer or else handwritten into the control chart. For the measurement the cell is placed in the compartment for the round cell and identified by the photo-meter via the bar code, and the measured absorbance is compared with the desired value. The absorbance is shown on the display and can be entered into the corresponding control chart.

The measurement of four cells for a given wavelength tests – in addition to the wavelength accuracy – also the linearity of the absorbance over the effective range.

The verification of the instrument, as it is required by DIN/ISO 9000 or GLP, can be easily performed by using the Spectroquant® PhotoCheck. The PhotoCheck hence offering the possibility to check the instrument. All of the corresponding documentation, required by these certification guidelines, is done by the photometer automatically.

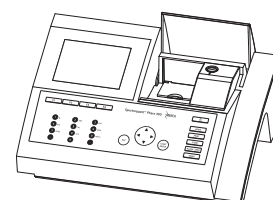

5.2.2 Checking the Overall System

Test for the overall system includes checking the working equipment and checking the handling operations.

The **overall system** can be checked using standard solutions of a known content, preferably with the Spectroquant® CombiCheck; this corresponds with the **AQA 2 mode** in the photometer.

Spectroquant® CombiCheck are ready-to-use standard solutions that in terms of the analyte concentration are finely adjusted to the individual test kits. They contain a mixture of several analytes that do not interfere with each other. The standard solution (R-1) is used in the same way as a sample. A double determination is recommended as a measure to diagnose any random errors.

The desired values with the permissible tolerances are already electronically stored with the method in the photometer.

In addition to the CombiCheck, it is also possible to use **CertiPUR® standard solutions** for this checking procedure. These contain 1000 mg of the respective analyte per liter of solution.

They can be diluted to different final concentrations, which should preferably lie approximately in the middle of the measuring range of the respective test kit. The table presented in Appendix 2 provides an overview of the available CombiCheck and ready-to-use standard solutions.

5 Analytical Quality Assurance (AQA)

Due to limited shelf-life characteristics, there are no CombiCheck or ready-to-use standard solutions for certain parameters. Appendix 3 is a compilation of **standard working procedures** necessary to make your own solutions of a defined concentration. This allows the control of parameters where there are no simple to prepare solutions available.


If the test for the overall system shows that all requirements are fulfilled, the individual results are flagged as AQA2. If not, an error message is given and the individual components of the instrument have to be checked in detail.

5.2.3 Checking the Pipettes



The **Spectroquant® PipeCheck** is used to check the pipettes. The pack contains cells filled with colour-dye concentrates. After the addition of a predefined volume of water using the pipette in question, the cell is measured against a corresponding reference cell also contained in the pack. The difference in the absorbance values of the measurement cell and reference cell may not exceed the tolerances given in the package insert. If the tolerances are exceeded, the instructions given in the section “Pipetting system” must be followed accordingly.

5.2.4 Checking Thermoreactors

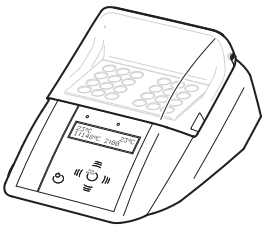

This is checked by means of the thermosensor. The thermoreactor is preheated as described in the Instructions for use. When the control lamp goes out, the temperature is measured in any one of the bores of the thermoreactor. The following desired temperatures must be achieved:

Block temperature 100 °C = desired temperature 100 ±3 °C
Block temperature 120 °C = desired temperature 120 ±3 °C
Block temperature 148 °C = desired temperature 148 ±3 °C

The even distribution of the temperature over all bores can also be documented using the thermosensor.

5 Analytical Quality Assurance (AQA)

5.2.5 Testing for Handling Errors

The user's own mode of operation must also be subjected to an exact analysis.

The following questions may serve as a guide in this regard:

- Is the test kit optimal for the measurement assignment in question?
- Is the test kit's measuring range suitable?
- Were the operating instructions for the test followed?
- Was the sample volume correct?
- Was the pipette handled properly?
- Was a new pipette tip used?
- Is the pH of the sample and measurement solution correct?
- Was the reaction time adhered to?
- Does the sample and reagent temperature lie within the correct range?
- Is the cell clean and free from scratches?
- Has the expiry date for the test kit been exceeded?

5.3 Determination of Sample Influences (matrix effects)

The influence of other substances contained in the sample may, under certain circumstances, be so great that their recovery rates lie in the region of several percent. It is recommended to check for any influence by using the addition solution contained in the Spectroquant® CombiCheck pack.

A defined quantity of the **addition solution** (R-2), which contains a known concentration of the respective analyte, is added to the sample and the recovery rate is determined. The following difference is then calculated:

$$\text{Result (sample + addition solution)} - \text{Result (sample)}$$

If the calculated difference is equal to the concentration of analyte of addition solution that was added, the recovery rate is 100 %. If the difference is less than 90 %, then a matrix interference is present.

5 Analytical Quality Assurance (AQA)

5.4 Definition of Errors

It is obvious that measurement results as a rule may be associated with errors. This applies equally to standardized methods of analysis (reference methods) and to routine analysis. The discovery and the minimization of errors must be the objective here.

A distinction is made between systematic errors and random errors.


Systematic errors are present when all the results of an analysis deviate from the true value with the same algebraic sign. Examples here include: a wrong sample volume, a wrong pH, a wrong reaction time, a sample-matrix influence, etc. Systematic errors thus affect the **accuracy** of the method of analysis.

Accuracy = Deviation of the measured concentration from the true concentration


Random errors manifest themselves in the form of a wide range of deviation of the results of a given sample. These can be kept to a minimum by ensuring good operating techniques and multiple determination with calculation of the mean values. Random errors make the result of the analysis unreliable; they influence the **precision**.

Precision = Dispersion of the results among each other


The following diagram illustrates the aspects of accuracy and precision:

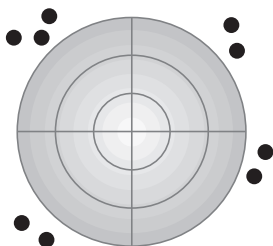


Accuracy: poor
Precision: poor
Major errors have been made!

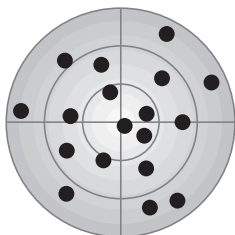


Accuracy: good
Precision: poor
Calculation of the mean values from at least three – or better even more – parallel determinations yields an approximation of the true value.

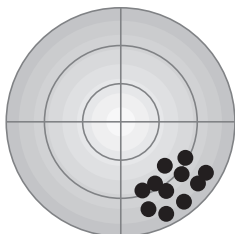

Accuracy: poor
Precision: good
The high degree of precision mistakenly indicates a correct value!

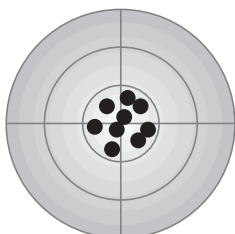

Accuracy: good
Precision: good
The ideal objective!

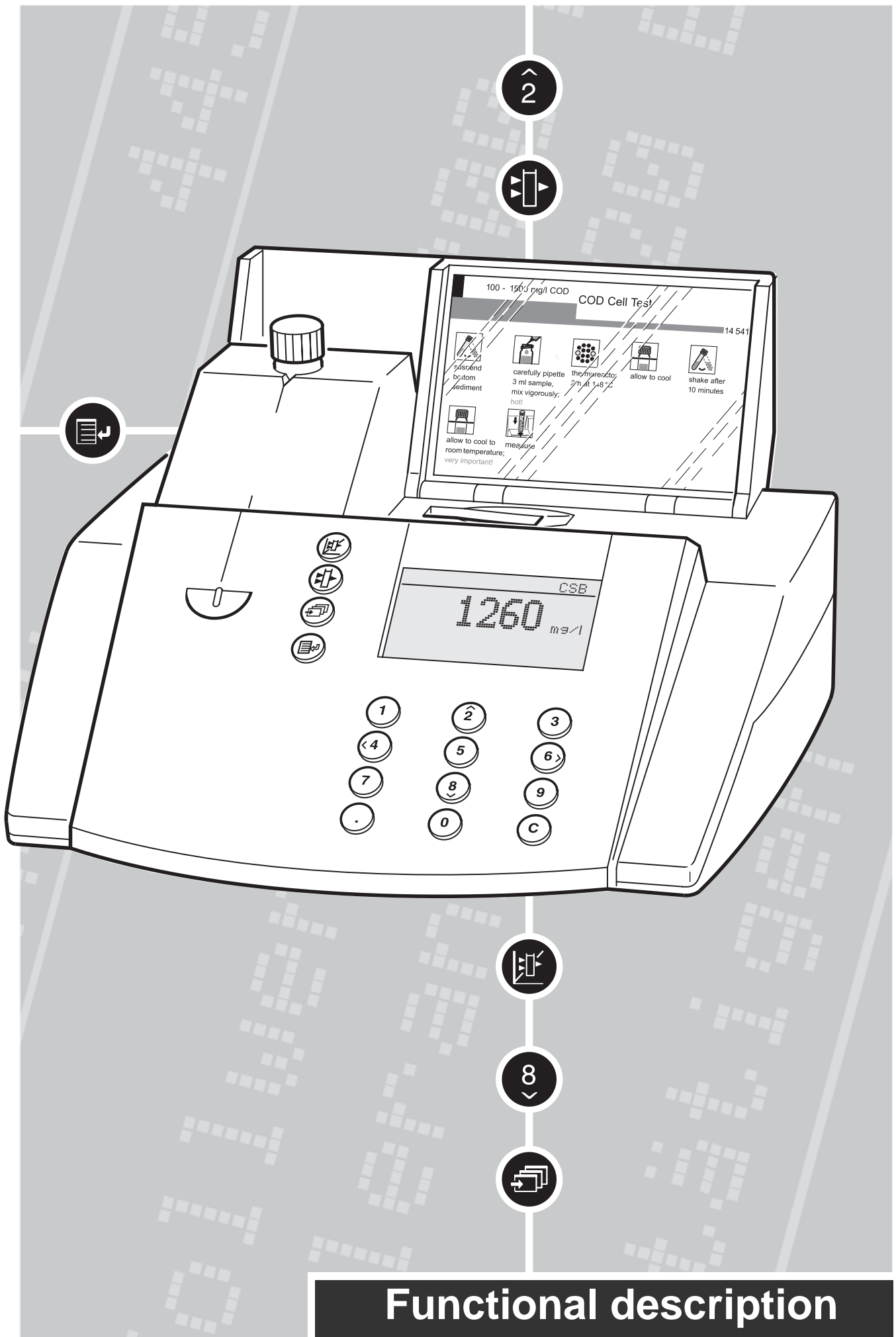

Functional description

General instructions

Notes on this operating manual

To ensure that you become rapidly acquainted with your photometer, the first chapter contains an overview and a short manual of the meter. The second chapter contains notes for the safe operation of the photometer.

Chapter 3 describes the commissioning of the photometer. The remaining chapters provide a comprehensive description of the functions and technical data of the photometer.

Symbols used



indicates notes that you must read – for your own safety, the safety of others and to protect your meter from being damaged.



indicates notes that draw your attention to special features.

Scope of delivery

- Photometer
- Power pack
- AQA MemoChip
- Product documentation

Warranty

The designated meter is covered by a warranty of 2 years from the date of purchase. The meter warranty extends to manufacturing faults that are determined within the period of warranty. The warranty excludes components that are replaced during maintenance, such as batteries, accumulators, lamps etc.

The warranty claim extends to restoring the meter to readiness for use but not, however, to any further claim for damages. Improper handling or unauthorized opening of the instrument invalidates any warranty claim.

To ascertain the warranty liability, return the meter and proof of purchase together with the date of purchase freight paid or prepaid.

Accuracy when going to press

The use of advanced technology and the high quality standard of our instruments are the result of continuous development. This may result in differences between this operating manual and your meter. We cannot guarantee that there are absolutely no errors in this manual. We are sure you will understand that we cannot accept any legal claims resulting from the data, figures or descriptions. The information in this manual is subject to change without notice.

1. Overview	5
1.1 Description of the operating elements	5
1.2 Identifying the connectors	5
1.3 Short manual	6
1.3.1 Measuring the concentration	6
1.3.2 Activating AQA2	6
1.4 Selecting and calling up the menu items	7
2. Safety	8
2.1 Authorized use	8
2.2 General instructions	8
2.2.1 Labeling of notes	8
2.2.2 Dangers of disregarding the safety instructions	8
2.2.3 Qualification of the personnel	8
2.2.4 Technical state of the meter	8
3. Commissioning	9
3.1 Preparing the photometer	9
3.2 Switching on the photometer	9
4. Measuring the Concentration	10
4.1 Measuring using cell tests	10
4.2 Measuring using reagent tests	11
4.3 Measuring using tests without barcode (manual method selection)	12
5. Measuring the Absorbance/Transmission	13
5.1 Switching to the Absorbance/ Transmission measuring mode	13
5.2 Measuring the absorbance or transmission	13
5.3 Measuring using cell tests	14
5.4 Measuring using reagent tests	15
5.5 Measuring using tests without barcode	15
6. Documentation	16
6.1 Resetting the number of the measured value	17
6.2 Activating the ID number	17
6.3 Download memory	18
6.4 Download of the methods list	21
7. Method Parameters	22
7.1 Citation form	23
7.1.1 Changing the citation form	23
7.1.2 Performing a difference measurement	24
7.2 Selecting the unit	26
7.3 Entering the dilution	27
8. Analytical Quality Assurance (AQA)	28
8.1 Activating AQA	28

Contents


8.1.1	Activating AQA using the AQA MemoChip	28
8.1.2	Activating AQA via the menu guide	28
8.1.3	Changing AQA intervals	30
8.1.4	Locking the system	32
8.1.5	Changing the password	33
8.1.6	Performing an AQA reset	33
8.2	Photometer monitoring (AQA1)	34
8.2.1	Entering PhotoCheck standards	34
8.2.2	Download of PhotoCheck standards	36
8.2.3	Erasing PhotoCheck standards	36
8.2.4	Performing Photometer monitoring	37
8.3	Total system monitoring with standard solutions (AQA2)	39
8.3.1	Entering standards	39
8.3.2	Output of standards	41
8.3.3	Erasing standards	42
8.3.4	Monitoring of the total system using standard solutions	42
9.	Kinetics	44
10.	Correction functions	45
10.1	Blank value	45
10.1.1	Activating the blank value measurement	46
10.1.2	Measuring the blank value	47
10.1.3	Erasing blank values	47
10.1.4	Recalling blank values	48
10.2	Reference absorbance	49
10.3	Turbidity correction	51
11.	Zero adjustment	52
12.	User-defined methods	54
12.1	Entering characteristics via the keyboard	55
12.2	Entering characteristics via the PC	57
12.3	Printing characteristics	58
12.4	Erasing user-defined methods	58
13.	Meter Setup	60
13.1	Selecting the language	60
13.2	Setting the date/time	61
13.3	Reset	61
13.4	System info	62
14.	Updating method data	63
15.	RS 232 C interface	65
15.1	Principle course of the remote control	65
15.2	Command list	65
15.3	Output format of measured values	66
15.4	Data transmission	66

Contents

15.5 Pin assignment	66
16. Maintenance	67
16.1 Changing the lamp	67
16.2 Actions to take if a cell is broken	67
17. Technical Data	68
18. What to do if...	71


1. Overview

1.1 Description of the operating elements

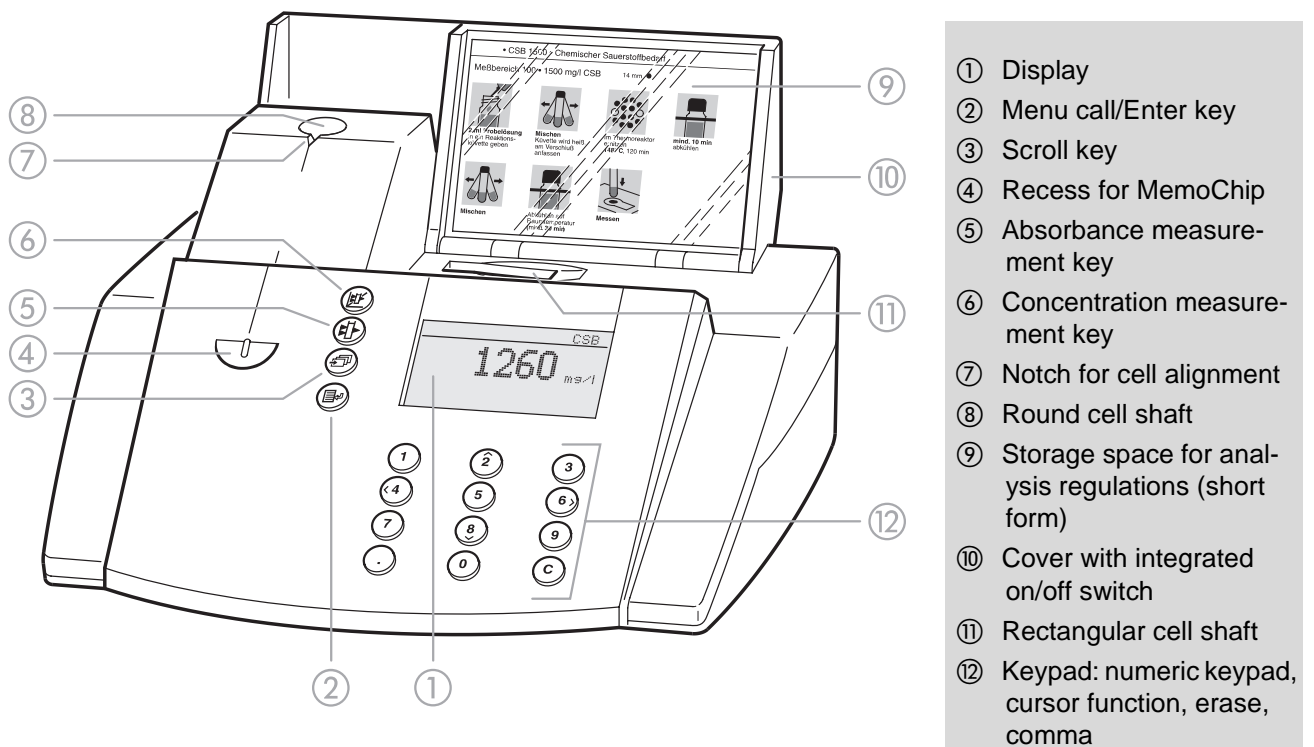

1.2 Identifying the connectors

- ① Power pack
- ② Connection for power pack
- ③ RS 232 interface

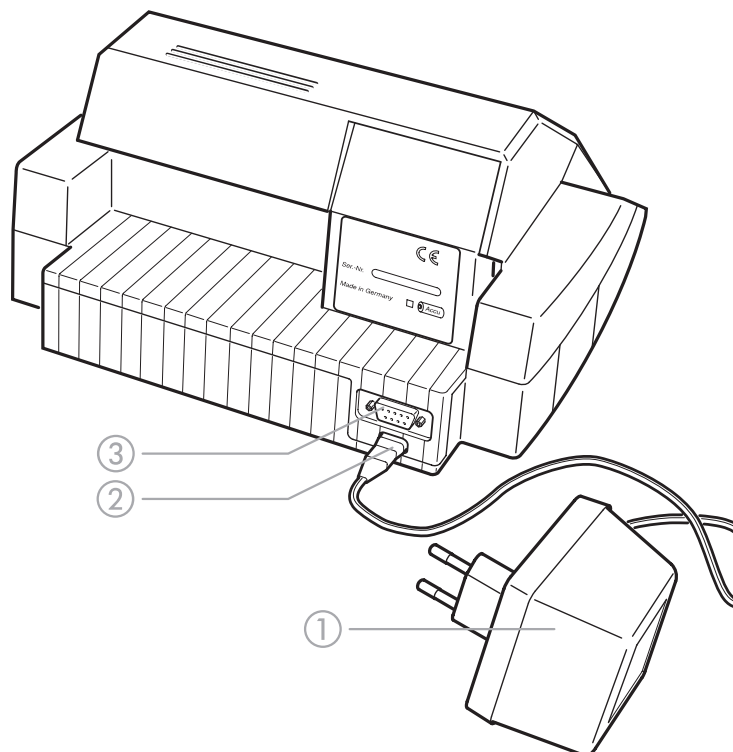


1. Overview

1.3 Short manual


The short manual lists all of the steps necessary to determine the concentration of a sample and to activate AQA2 at a glance.

1.3.1 Measuring the concentration


– To switch on the photometer, open the cover.
The photometer performs a check (*Self-Check*) of the entire system and then switches automatically to the *concentration measuring mode*.

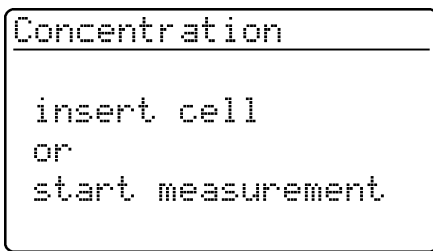

Measuring mode, *concentration*

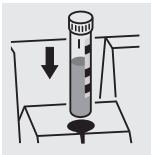


or

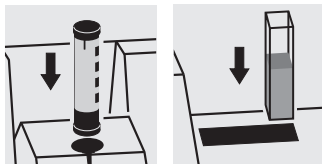

- Insert the round cell with barcode in the round cell shaft until it clicks into place, or
- Insert the AutoSelector in the round cell shaft and the rectangular cell in the rectangular cell shaft. Align the line mark to the notch of the photometer. The message *measuring... appears*.



If the *select method* menu is displayed, align the line mark of the round cell or AutoSelector to the notch of the photometer.

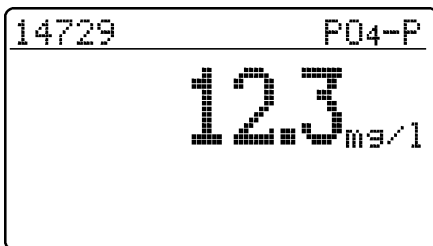


The measured value appears on the display. Measured values outside the specified measuring range are output in small numerals. Repeat the measurement:

- Press .

1.3.2 Activating AQA2

– Place the AQA MemoChip in the recess on the photometer.
The following display appears:

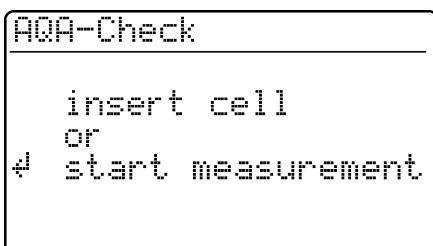


The action of placing the AQA MemoChip in the recess directly activates the AQA check without having to press any key.



Standard concentrations and tolerances are listed in the table "Spectroquant® CombiCheck and standard solutions" in the part, "General information".

1. Overview

1.4 Selecting and calling up the menu items

- To switch on the photometer, open the cover.
- Press .

The following display appears:

```
Setup
-----
▶documentation
  method Parameter
  Kinetic
  Meter Setup
```

The following display appears:

```
Setup
-----
  documentation
  method Parameter
  Kinetic
  ▶Meter Setup
```

```
meter setup
-----
  return
  ▶AQA Functions
  Correction Funct.
  adjust zero
  user def. methods
```


Example:

The *documentation* menu item is preselected in the *setup* menu (▶).





Select a menu item, e. g. *meter setup*:

- Press  or  .






The *meter setup* menu item is preselected (▶).

- Call up the *meter setup* submenu by pressing .



The required menu item is

- selected using  or  .
- called up using .

Selection lists:

- Changes to the settings are accepted after confirmation by pressing .
- Current settings are marked by "◆".
- Change to other configuration levels by
 - Selecting the menu item, *return*
 - Pressing 
- Scroll with  or  .

Character input:

- via the numeric keypad or by using ,
- the character to be input is shown in reverse video.
- Confirm each input with .

2. Safety

This operating manual contains basic instructions to be followed in the commissioning, operation and maintenance of the meter. Consequently, all responsible personnel must read this operating manual before

working with the meter. The operating manual must always be available in the vicinity of the meter.

2.1 Authorized use

The photometer is authorized exclusively for analyzing substances in water and aqueous solutions using round cells or rectangular cells (special optical glass).

Observe the technical specifications of the cells according to chapter 17 TECHNICAL DATA. Any other use is considered **unauthorized**.

2.2 General instructions

The photometer is constructed and tested according to the EN 61010-1 safety regulations for electronic measuring instruments. It left the factory in a safe and secure technical condition.

The smooth functioning and operational safety of the photometer can only be guaranteed under the climatic conditions specified in chapter 17 TECHNICAL DATA of this operating manual.

Opening the photometer or adjustment, maintenance and repair work must only be performed by personnel authorized by the manufacturer.

The only exceptions to this are the activities described in chapter 16 MAINTENANCE. Non-compliance results

in the loss of warranty claims.

Follow the points listed below when operating the photometer:

- Follow local safety and accident prevention regulations.
- Observe the enclosed instructions concerning reagents and accessories.
- Observe the regulations when dealing with dangerous substances.
- Follow the operating instructions at the workplace.
- Use only original spare parts.

2.2.1 Labeling of notes



indicates notes that you must read – for your own safety, the safety of others and to protect your meter from being damaged.



indicates notes that draw your attention to special features.

2.2.2 Dangers of disregarding the safety instructions

Disregarding the safety instructions can adversely affect the safety of both the user and the environment as well as the equipment.

Non-compliance with the safety instructions will result in the loss of any warranty claims.

2.2.3 Qualification of the personnel

The personnel responsible for the commissioning, operation and maintenance must have the necessary qualifications for this work. If the personnel do not have the required skills they have to be instructed.

Furthermore, it must be ensured that the personnel read and completely understand the present operating manual.

2.2.4 Technical state of the meter

It is the responsibility of the operator to continuously observe the overall technical condition (externally recognizable deficits and damage as well as alterations to the operational behavior) of the meter.

If safe operation is no longer possible, the equipment must be taken out of service and secured against inadvertent operation.

Safe operation is no longer possible if

- the equipment has been damaged in transport
- the equipment has been stored under adverse conditions for a lengthy period of time
- the equipment is visibly damaged
- the equipment no longer operates as prescribed.

If you are in any doubt, please contact the supplier of the photometer.

3. Commissioning

The photometer operates at an environmental temperature of +5 °C to +40 °C. During transport from cold to warm surroundings, condensation can form resulting in the malfunction of the meter.

Before putting the photometer into service, wait until it has adapted to the new environmental conditions (see also chapter 17 TECHNICAL DATA).

3.1 Preparing the photometer

- Place the photometer on a hard, flat surface and protect it against intensive light and heat.

Line operation


- Plug the original power pack into the socket on the photometer
- Plug the power pack into the line socket
- Switch on the photometer (open the cover).

Battery operation

- Charge the battery for approx. 5 hours before the initial commissioning. To do this:
 - Plug the original power pack into the socket on the photometer
 - Plug the power pack into the line socket and then the battery will be charged.


During battery operation or when the meter is at a standstill for longer periods of time, the battery runs down. This can result in your photometer no longer being ready for operation.

When the following symbol is displayed, charge the

battery: 


3.2 Switching on the photometer

- To switch on the photometer, open the cover. The photometer performs a check (*Self-Check*) of the entire system and then switches automatically to the *concentration* measuring mode.

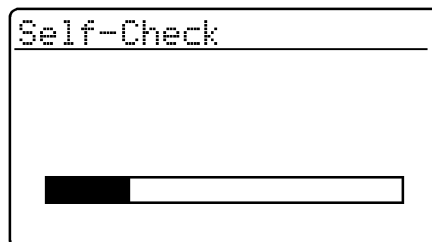


Self-check of the photometer


After approx.
5 s:

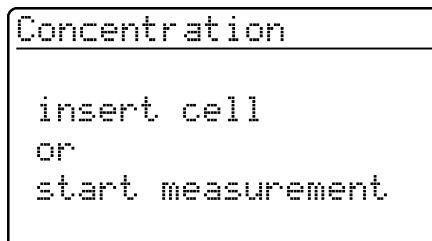

Automatic change to the measuring mode,
concentration

4. Measuring the Concentration



- Call up the *concentration* measuring mode by actuating .

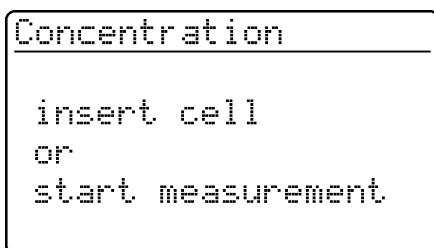

```
Concentration
-----
insert cell
or
start measurement
```

Measuring mode, *concentration*

4.1 Measuring using cell tests

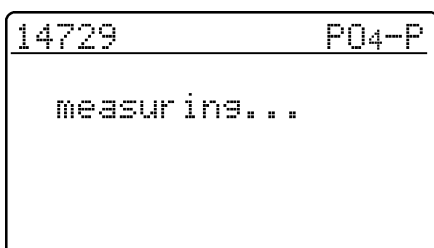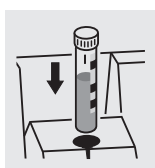


```
14729 P04-P
-----
measuring...
```

- Insert the round cell with barcode into the round cell shaft until it clicks into place.
Align the line mark to the notch of the photometer.

The photometer reads the barcode of the round cell and automatically selects the relevant method.

After approx.
2 s:

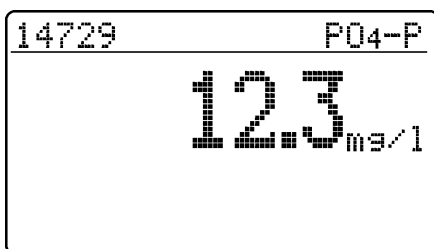

```
14729 P04-P
-----
12.3 mg/l
```


The measured value appears on the display.



If the *select method* menu is displayed, align the line mark of the round cell or AutoSelector to the notch of the photometer.


4. Measuring the Concentration

4.2 Measuring using reagent tests

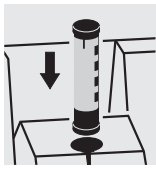

14773 NO₃-N
insert cell
or
start measurement

- Insert the AutoSelector into the round cell shaft. Align the line mark to the notch of the photometer. The photometer reads the barcode and automatically selects the relevant method.

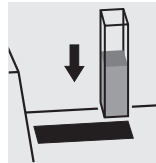





14773 NO₃-N
measuring...

- Insert the rectangular cell (vertically, to the right of the slot). The measurement starts automatically. Observe the position mark at the cell shaft.



If the *select method* menu is displayed, insert the corresponding AutoSelector into the round cell shaft.

If the "cell type" query appears, select the relevant rectangular cell (10 mm, 20 mm, 50 mm) with  or  , confirm with .

After approx.
2 s:

14773 NO₃-N
12.3 mg/l



The measured value appears on the display.

4. Measuring the Concentration


4.3 Measuring using tests without barcode (manual method selection)


When measuring using cell tests without barcode or reagent tests without AutoSelector, the method must be selected manually.

The sequence of operations also applies to user-defined methods.

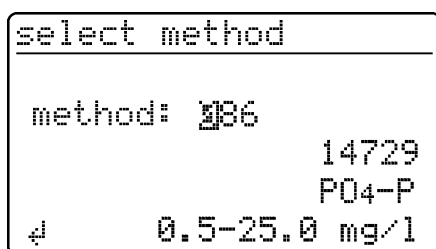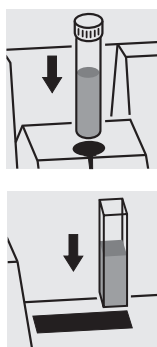


The last method set up manually appears on the display.

- Enter the required method number via the numeric keypad
- Confirm with .

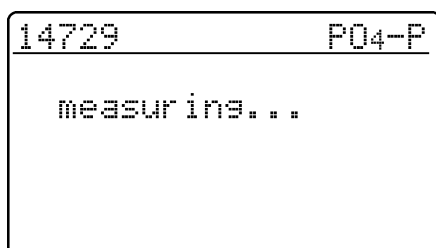


After approx.
2 s:

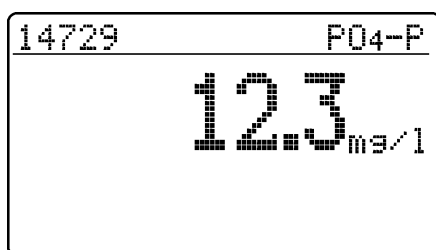

The measured value appears on the display.

5. Measuring the Absorbance/Transmission

5.1 Switching to the Absorbance/ Transmission measuring mode

- Call up the *setup* measuring mode by actuating .

```
Setup
-----
documentation
method Parameter
Kinetic
▸abs./trm. %
Meter Setup
```


- In the *setup* menu, call up the *abs./trm. %* submenu.

```
▸absorbance *
transmission
return
```

Selection of the measuring mode:

- *absorbance*
- *transmission*

5.2 Measuring the absorbance or transmission

- Call up the *absorbance* or *transmission* measuring mode (depending on the selection in the *abs./trm. %* menu) by actuating .

```
Absorbance
-----
insert cell
or
start measurement
```

Measuring mode, *absorbance*

```
transmission
-----
insert cell
or
start measurement
```



Measuring mode, *transmission*

i The transmission measurement is not described separately in the following example as it operates in exactly the same way as the absorbance measurement. However, the result of the measurement is displayed as % Transmission instead of A for Absorbance.


i A measured reference absorbance is also effective in the measuring mode, *transmission*. It is displayed as reference absorbance.

5. Measuring the Absorbance/Transmission

5.3 Measuring using cell tests

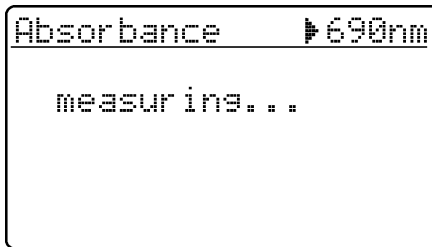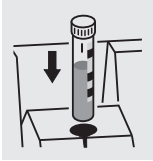



- Insert the round cell with barcode into the round cell shaft until it clicks into place.
Align the line mark to the notch of the photometer.

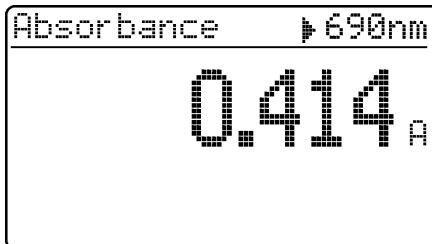


- The measured value for the wavelength displayed at the top right appears. This measured value is automatically stored.
If necessary, call up further wavelengths:

- with  or  .

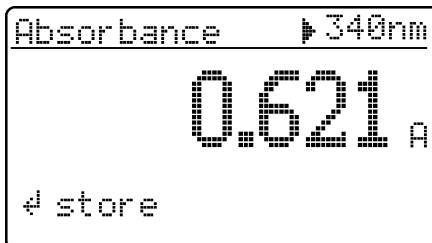



- The measured value for the selected wavelength appears and can be stored and output to the interface with .

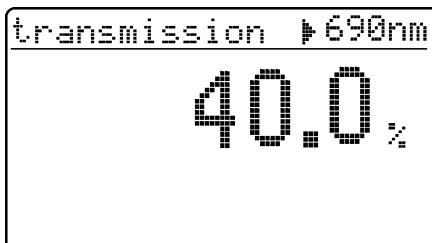



- Sample display for transmission measurement

5. Measuring the Absorbance/Transmission

5.4 Measuring using reagent tests

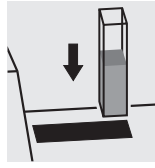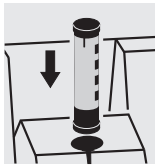


- Insert the AutoSelector in the round cell shaft and the rectangular cell in the rectangular cell shaft. Align the line mark to the notch of the photometer.

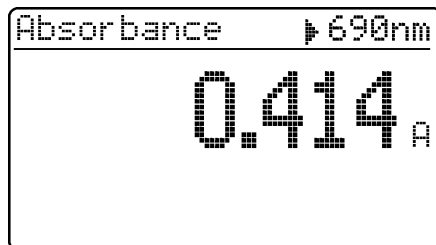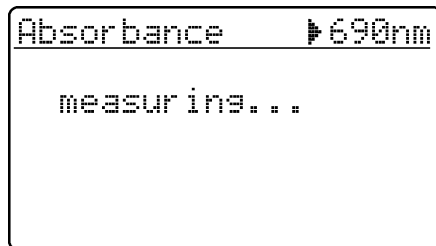

The measured value for the wavelength displayed at the top right appears. This measured value is automatically stored.



If necessary, call up further wavelengths:

- with  or  .

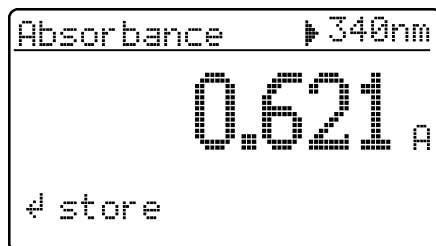

The measured value for the selected wavelength appears and can be stored.

5.5 Measuring using tests without barcode

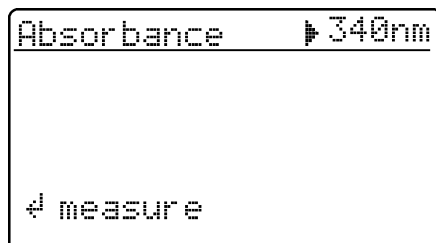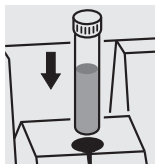

The last wavelength measured appears on the display.

- Select the wavelength:


with  or  .

- Start the measurement:

Press .

6. Documentation

The measured values can be documented as follows:

- Storage in the measured value memory
 - Output to a connected printer via the serial interface (automatic when a printer is connected)
 - Transmission to a PC for further processing (by using the relevant software, e.g. Multi/ACHATII or – less conveniently – by means of a terminal program).
- To switch on the photometer, open the cover.
- Press .

The following display appears:

```
Setup
-----
▶documentation
  method parameter
  Kinetic
  Meter Setup
```

```
documentation
-----
▶no. of meas. value
  I.D. number
  download memory
  output methods
  return
```

- Call up the *documentation* menu with .

The following functions can be selected:

- *no. of meas. value*
 - reset the number
- *I. D. number*
 - On / Off
- *download memory*
 - total
 - from date
 - with I. D. no. XXX
 - for method XXX
 - AQA
- *output methods*
 - all
 - user def. methods

The current settings are marked by "◆" in the selection lists of the respective submenus.

6. Documentation





6.1 Resetting the number of the measured value

```
documentation
-----
▶no. of meas. value
  I.D. number
  download memory
  output methods
  return
```

```
no. of meas. value
-----
reset number:
▶Yes  +
  No
  return
```

- Call up the *no. of meas. value* submenu.

- *yes*
The numbering of the measured values starts again with 001 (default)
- *no*
Consecutive numbering of the measured values (from 001 to 999)

- Select the menu item with  or  
- Confirm with .

6.2 Activating the ID number

If the *I. D. number* function is active, a sequence of up to 6 alphanumeric characters (I. D. number) is allocated to a measurement (e. g. sample location, customer number, river kilometer).



```
documentation
-----
no. of meas. value
▶I.D. number
  download memory
  output methods
  return
```

```
I.D. number
-----
▶Off  +
  on
  return
```

```
I.D. number
-----
Off
▶on  +
  return
```


- Call up the *I. D. number* submenu.

- *off*
It is not possible to enter an ID number (default).
- *on*
Input an ID number for each measured value.

- Select the menu item, *on*, with 
- Confirm with .


6. Documentation

Measuring using the activated "ID number" function

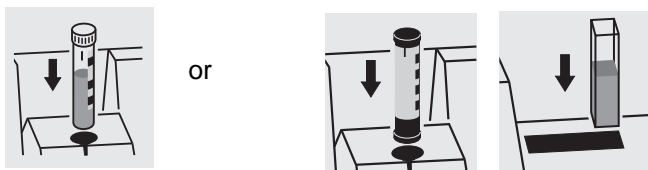

- Call up the *concentration* measuring mode
- Insert and align the round cell, or
- Insert the AutoSelector and rectangular cell.

The following display appears:

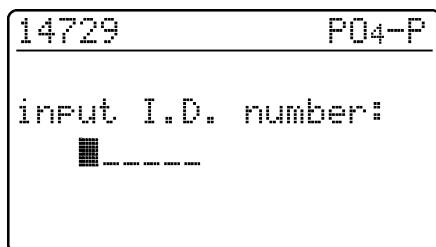

- Enter the required ID number
Default:
The ID number that was last entered
(initially, there are underscores).

The input is made as follows:



- Capital letters with
- Numerals via the numeric keypad
- Confirm each time with . Confirmation without entering a character results in a space.

To erase incorrect inputs:

- By pressing the last input is erased.

6.3 Download memory

The measured value storage can be selectively downloaded to either the display or serial interface. The selection of the output medium is made after the specification of the sorting criteria.

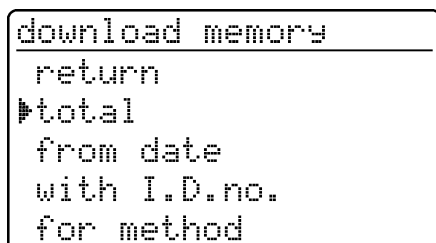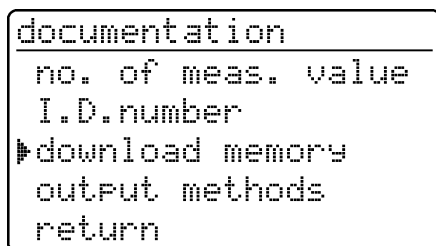

- Call up the *download memory* submenu.

The *download memory* menu item only appears after at least one measurement has been performed.

The following sorting criteria can be set:

- *total* – all stored measured values
- *from date* – all measured values from a special date
- *with I. D. no.* – all measured values of a specific ID number
- *for method* – all measured values of a special method
- *AQA* – all measurements of a specific method monitored with AQA.





- Select the menu item with
- Confirm with .

6. Documentation

Selecting "total"



```
download memory
└to display
  to Printer/PC
  return
```

Select the output medium:

- *to display*
- *to printer/PC* (serial interface).
- Select the menu item with  or  
- Confirm with  to start the memory download.





Selecting "from date"

```
download memory
from date: 20.02.98
↵
```

- Input the date via the numeric keypad
- Erase the input using 
- Confirm with .

```
download memory
└to display
  to Printer/PC
  return
```


Select the output medium:

- *to display*
- *to printer/PC* (serial interface).
- Select the menu item with  or  
- Confirm with  to start the memory download.





6. Documentation

Selecting "with I. D. no."

```
download memory
-----
input I.D.number:
FEED
```


- Enter the ID number
Default:
The last ID number entered is preselected.
- Confirm each time with .

```
download memory
-----
to display
to printer/PC
return
```





- Select the output medium:
- *to display*
 - *to printer/PC* (serial interface).
- Select the menu item with  or  .
 - Confirm with  to start the memory download.

Selecting "for method"

```
download memory
-----
method: 086
                14729
                P04-P
↵          0.5-25.0 mg/l
```


- Input the method
Default:
The last method selected.
- Confirm with  to start the memory download.

```
download memory
-----
to display
to printer/PC
return
```

- Select the output medium:
- *to display*
 - *to printer/PC* (serial interface).
- Select the menu item with  or  .
 - Confirm with  to start the memory download.

Selecting "AQA"

```
download memory
-----
method: 086
                14729
                P04-P
↵          0.5-25.0 mg/l
```

- Input the method
Default:
The last method selected.
- Confirm with  to start the memory download.

6. Documentation

Memory download to display

```
download memory
-----
009  07.05.97 17:24
Feed      14554 Ni
          3.66 mg/l
↓ return  AQA2
```

Each data record appears individually on the display beginning with the data record just measured.

The display shows:

- no. of meas. value
- date/time
- I. D. number
- method designation
- citation
- meas. value
- unit
- Where necessary, AQA ID, e.g. AQA2.


– Scroll with  or  .

Memory download to printer/PC

```
download memory
data transmission
runs:
          121
↓ cancel
```

Memory download to the serial interface:

- Display of the transmitted no. of measured value (continuation display) beginning with the last measured value.

– Cancel with .

Sample printout:

003	14541	10.02.98	11:56:33	t	80	mg/l	COD
002	14541	10.02.98	11:54:21	t	70	mg/l	COD
001	14729	03.02.98	18:30:53	*	0.3	mg/l	PO4-P

6.4 Download of the methods list

The stored methods are downloaded to the printer/PC via the serial interface.

```
documentation
-----
no. of meas.value
I.D. number
download memory
▶output methods
return
```

– Call up the *output methods* submenu.

```
output methods
-----
▶all
user def.
return
```

The following parameters can be set:

- *all* – Download of all stored methods
- *user def.* – Download of user-defined methods.



The menu item, *user def.*, appears only after user-defined methods have been input.

– Select the menu item with  or  .

– Start the download with .

7. Method Parameters

The following parameters can be set in the *method parameters* menu:

- citation
- unit
- dilution

– To switch on the photometer, open the cover.

– Press .

The following display appears:


```
Setup
-----
documentation
▶method Parameter
Kinetic
Meter Setup
```

```
method parameter
-----
method: 086
                14729
                PO4-P
↵          0.5-25.0 mg/l
```

```
method parameter
-----
▶Citation
unit
Dilution
return
```

– Call up the *method parameters* submenu.

– Input the method number

– Confirm with .

– Select the menu item with  or  .

– Call up the parameter by pressing .

7. Method Parameters

7.1 Citation form

7.1.1 Changing the citation form

Example:

Change the citation form from "NH₄-N" to "NH₄".

```
method parameter
└─ Citation
  unit
  Dilution
  return
```


– Call up the *citation* submenu.

```
Citation 14739
└─ NH4-N ⚡
  NH4
  return
```

The current setting: NH₄-N (⚡).

```
Citation 14739
  NH4-N ⚡
└─ NH4
  return
```

– Using , scroll to NH₄

– Confirm with .


```
Citation 14739
  NH4-N
└─ NH4 ⚡
  return
```

– Citation form NH₄ is set (⚡).

7. Method Parameters

7.1.2 Performing a difference measurement


Difference measuring is possible for some methods (e.g. Iron II/III, Ca-/Mg Hardness).

 For more information on this, see part, "Analysis specifications".

Example:

Determination of iron (II) and iron (III).

```
method parameter
-----
method: 106
                        14896
                        Fe
↓      1.0-50.0 mg/l
```



- Enter method 106
- Confirm with .

```
method parameter
-----
▶Citation
unit
Dilution
return
```

- Call up the *citation* menu item.


```
Citation      14896
-----
▶Fe  +
FeII,FeIII Δ
return
```


The current setting: Fe

- Using  scroll to Fe II, Fe III Δ
- Confirm with .

```
Citation      14896
-----
Fe
▶FeII,FeIII Δ +
return
```

- Citation form Fe II, Fe III Δ (+) is set.

- Change to measuring by pressing .

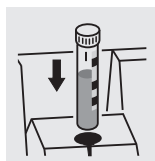

```
14896      FeII,FeIII Δ
-----
Σ Fe
measuring...
```

- Start the 1st measurement by inserting cell 1.


7. Method Parameters


After approx.
2 s:

14896	FeII,FeIII Δ
3.2 mg/l	
Σ Fe	
\leftarrow FeII	

The 1st measured value appears on the display: Σ Fe.

– Remove cell 1

– Press .

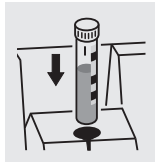

14896	FeII,FeIII Δ
FeII	
measuring...	

– Start the 2nd measurement by inserting cell 2.

After approx.
2 s:

14896	FeII,FeIII Δ
2.1 mg/l	
FeII	
\leftarrow FeII,FeIII	

The 2nd measured value appears on the display:
Iron II.

– Continue to the display of both measured values

using .

14896	FeII,FeIII Δ
FeII	2.1 mg/l
FeIII	1.1 mg/l

Display of both measured values as a summary.

7. Method Parameters

7.2 Selecting the unit



The preset unit is "mg/l".
It can be changed to "mmol/l".

```
method parameter
-----
Citation
unit
Dilution
return
```

- Call up the *unit* submenu.

```
unit 14729
-----
mg/l  +
mmol/l
return
```

The current setting: *mg/l* (⊕)

- Using  scroll to *mmol/l*
- Confirm with .

```
unit 14729
-----
mg/l  +
mmol/l
return
```

- Unit *mmol/l* (⊕) is set.

7. Method Parameters

7.3 Entering the dilution

Diluting a sample with **distilled water** enables the measuring range to be extended.

The photometer automatically includes any dilution factor entered when calculating the measured value.

Values from 0 to 99 can be entered as dilution factors.

Example:


Enter dilution 1:20 (i.e. 1 part sample + 19 parts distilled water).


```
method parameter
-----
Citation
unit
▶Dilution
return
```

- Call up the *dilution* submenu.



```
Dilution 14729
-----
sample + water:
  1 + 00
↵
```

The current setting: 1:00.

- Input the dilution factor 19 via the numeric keypad
- Confirm with .

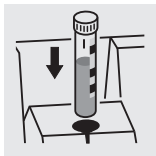

or

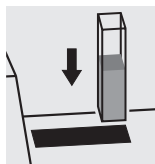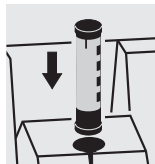

- Call up the *concentration* measuring mode
- Insert the cell
- Confirm the method.

```
14729 PO4-P
-----
 127 mg/l
      U:1+19
```

- The measured value is displayed together with the dilution factor.

Repeat measurements:

- Insert the cell, or start the measurement
- Confirm the dilution factor.

Erase the dilution factor by:

- changing the method
- switching off the photometer
- entering 00 as dilution number.

8. Analytical Quality Assurance (AQA)

Analytical quality assurance (AQA) can be performed in two steps:

- **AQA1** – Photometer monitoring
- **AQA2** – Total system monitoring with standard solutions



The total system monitoring (AQA2) is a method-specific check using standard solutions.

If this is performed successfully, it also includes photometer monitoring (AQA1).

See also part "General information" for further information on Analytical Quality Assurance (AQA).

The AQA mode must be activated in the photometer. In the delivery state it is switched off.

The AQA mode is activated:

- by inserting the AQA MemoChip
 - monitoring of the total system using standard solutions (AQA2)
- by using a menu to select
 - monitoring of the photometer (AQA1)
 - monitoring of the total system using standard solutions (AQA2)

8.1 Activating AQA

- To switch on the photometer, open the cover.

8.1.1 Activating AQA using the AQA MemoChip

- Place the AQA MemoChip in the recess on the photometer.

The following display appears:



```
AQA-Check
-----
insert cell
or
↵ start measurement
```



The action of placing the AQA MemoChip in the recess directly activates the AQA2 check (see section 8.3.4).

8.1.2 Activating AQA via the menu guide

- Press

```
Setup
-----
documentation
method parameter
Kinetic
▶Meter Setup
```

- Call up the *meter setup* submenu.

```
meter setup
-----
return
▶AQA Functions
Correction Funct.
adjust zero
user def. methods
```

The *meter setup* submenu appears with the *AQA functions* menu item preselected.

- Confirm with


8. Analytical Quality Assurance (AQA)

A password request appears:

```
AQA-Configuration
-----
input password:

  0 0 0 0
```

A separate password protects settings of the AQA-configuration against unauthorized access (Changing the password see section 8.1.5).

- Input the password:
Only **numeric** characters are allowed.
Default: 0000
- Confirm with .

If the input was incorrect:

```
AQA-Configuration
-----
wrong password
```

- Repeat the input.



If you have forgotten the password, contact the service department.

After the password has been successfully input, the *AQA configuration* submenu appears:


```
AQA-Configuration
-----
return
▶AQA-Mode
  AQA-Standards
  AQA-Intervals
  System locked
```

- Call up the *AQA mode* function.

```
AQA-Mode
-----
▶off *
  n weeks
  n measurements
  return
```

Default: *off* (no monitoring)

- Select *AQA mode*:
 - *off*
 - *n weeks*
 - *n measurements*

- Confirm with .

```
Setup
-----
documentation
method parameter
▶AQA-Check
  Kinetic
  Meter Setup
```

- In the *setup* menu, call up the *AQA check* submenu.

```
AQA-Check
-----
Meter
▶system
  return
```

Selection of the AQA mode:

- *meter*
- *system*



The menu item, *meter*, only appears after the corresponding PhotoCheck standards have been input (see section 8.2.1).

8. Analytical Quality Assurance (AQA)

8.1.3 Changing AQA intervals

AQA intervals specify the interval between two AQA checks. A fixed time interval (*n weeks*) or a number of measurements (*n measurements*) can be specified as the interval.

The respective values that were input remain stored even if they are not activated.

Additionally, two separate intervals can be set up for both photometer monitoring (AQA1) and system monitoring (AQA2).



For the total system monitoring (AQA2), a change of the time interval (*n weeks*) even retroactively applies to monitoring processes that are already running.

Changing the number of measurements (*n measurements*) does not affect monitoring processes already running.

Thus, individual numbers of measurements can be set for different methods.

```
AQA-Configuration
-----
return
AQA-Mode
AQA-Standards
▶AQA-Intervals
System locked
```

After an interval has expired, the following consequences become effective:

- Warning and loss of AQA identification
- Locking of the method for concentration measurements (as long as the locking is active).

Setting ranges:

- Photometer monitoring (AQA1):
 - 1 to 52 weeks (default: 12 weeks) or
 - 1 to 9999 measurements (default: 1500)
- Monitoring of the total system using standard solutions (AQA2):
 - 1 to 52 weeks (default: 4 weeks) or
 - 1 to 9999 measurements (default: 100)



With the *n measurements* setting, a difference measurement (see section 7.1.2) is counted as one measurement only.

- In the *AQA configuration* menu, call up the *AQA intervals* submenu.
According to the selection in the *AQA mode* menu, a fixed time interval (*n weeks*) or a number of measurements (*n measurements*) is set in the *AQA intervals* menu.



If the *AQA mode* function is switched off, the *AQA intervals* submenu is not available.

AQA interval, "n weeks"

The AQA interval, *n weeks*, is only effective if the *n weeks* setting is active for the *AQA mode* function.

The specified number of *n weeks* applies to:

- the photometer with AQA1
- all methods with AQA2.

- In the *AQA intervals* menu, call up the *n weeks* submenu.

```
AQA-Intervals
-----
AQA-Meter: 04 w
AQA-System: 04 w
← confirm
```

- To return without change, press three times
- Enter the time interval for *AQA meter* via the numeric keypad, confirm with .

8. Analytical Quality Assurance (AQA)

AQA interval, "n measurements"

The AQA interval, *n measurements*, is only effective if the *n measurements* setting is active for the *AQA mode* function.

The AQA2 check starts the monitoring for one method at a time.

The specified number of *measurements* applies to:

- the instrument with AQA1 (total number of measurements performed, independent of whether AQA2 is active for some parameters)
- each method an AQA check will then be performed for with AQA2.
Thus, it is possible to define individual numbers of measurements for different methods.
The measurements are counted separately for each monitored method.

The monitoring intervals of AQA2 monitoring processes already started for other methods are not affected by changing the number of *measurements*. Thus the number of *measurements* can be set for further methods no matter which monitoring processes were started before.






When an AQA2 check is performed, the number of *measurements* last set in the *AQA intervals* menu is automatically taken over.

Therefore, you should check and, if necessary, change the currently set number of *measurements* before each AQA2 check.

The currently set number of *measurements* for the AQA2 check is saved for the active method and output in the report individually (section 8.3.4).

- In the *AQA intervals* menu, call up the *n measurements* submenu.

```
AQA-Intervals
-----
AQA meter:
  1500 measurements
AQA system:
  0100 measurements
* confirm
```

- To return without change, press  three times
- Enter the number of measurements for *AQA meter* via the numeric keypad, confirm with .
- Enter the number of measurements for *AQA system* via the numeric keypad, confirm with .

8. Analytical Quality Assurance (AQA)

8.1.4 Locking the system

The function *system locked* is effective if, for a monitored method,

- no AQA check was performed,
- the AQA check “system” has expired.


As a result, a concentration measurement is not possible for this method.

```
AQA-Configuration
return
AQA-Mode
AQA-Standards
AQA-Intervals
▶System locked
```

```
System locked
off
▶on      +
return
```

– Call up the *system locked* submenu.

– Select the menu item with  or  .


– Confirm with .

8. Analytical Quality Assurance (AQA)


8.1.5 Changing the password

When delivered, the default password is 0000. This password can be changed as follows:


```
AQA-Configuration
-----
AQA-Standards
AQA-Intervals
System locked
▶change Password
reset
```

- Call up the *change password* submenu.
- Confirm with .

```
AQA-Password
-----
input Password:
(0000)
  0 0 0
```

- Input the required password, e.g. 0100, via the numeric keypad
- Confirm with .


```
AQA-Password
-----
confirm Password:
(0100)
  0 0 0
```

- Input the password once again:
- Confirm with .

8.1.6 Performing an AQA reset


If the Analytical Quality Assurance is to be switched off completely or reset to the delivery state, this can be made via the *reset* function in the *AQA configuration* submenu.

```
AQA-Configuration
-----
AQA-Intervals
System locked
change Password
▶reset
return
```

- Call up the *reset* submenu
- Confirm with .

```
AQA-Configuration
-----

▶reset
cancel
```

- Select the *reset* menu item
 - Confirm with .
- An AQA reset is performed.


8. Analytical Quality Assurance (AQA)

8.2 Photometer monitoring (AQA1)

8.2.1 Entering PhotoCheck standards



A Spectroquant® PhotoCheck is required to perform the photometer monitoring (AQA1). **At least 1 standard** must be input. We recommend, however, to input all available standards.

- Press  to call up the *setup* menu
- Call up the *meter setup* submenu.
- Call up the *AQA functions* submenu.
- Input the password
- Call up the *AQA standards* submenu and the following display appears:

```
AQA-Standards
├─PhotoCheck
│  standard solution
│  return
```

```
PhotoCheck-Standards
├─input
│  output
│  erase
│  return
```

- Call up the *PhotoCheck* submenu.

Select between

- *input*
Input the theoretical value (absorbance) from the lot certificate of Spectroquant® PhotoCheck
- *output*
Print/display theoretical values
- *erase*
Erase theoretical values.







The *erase* and *output* menu items only appear after at least one standard has been input.

8. Analytical Quality Assurance (AQA)


Example:

445-1 nm, theoretical value (absorbance) 0.200,
admissible tolerance ± 0.020

```
PhotoCheck-Standards
return
▶445-1
  445-2
  445-3
  445-4
```


- Select with  or  
- Quit via the menu item, *return*
- Confirm with .

```
PhotoCheck          445-1
theor.val.: 0.200 A
#confirm
```

- Input the theoretical value, 445-1
- Confirm with .


If the standard is already stored, this value appears on the display.

```
PhotoCheck          445-1
theor.val.: 0.200 A
Tolerance: ±0.020 A
#confirm
```

- Input the tolerance via the numeric keypad
- Confirm with .

```
PhotoCheck-Standards
return
▶445-1 ✓
  445-2
  445-3
  445-4
```

PhotoCheck standard 445-1 is input.

- Select the next one with 
- Input all PhotoCheck standards in this way.

8. Analytical Quality Assurance (AQA)

8.2.2 Download of PhotoCheck standards





```
PhotoCheck-Standards
input
output
erases
return
```

- In the *PhotoCheck standards* submenu, call up the *output* menu item.

```
download PhotoCheck
to display
to printer/PC
return
```

Select the output medium:

- to display
- to printer/PC (serial interface).

- Select with  or  
- Confirm with  to start the download.

Example: Report output

AQA check meter		AQA1		
26.08.97		13:19		
AQA interval		12 weeks		
test sol.	unit	theor. val.	tolerance	AQA date
445-1	A	0.200	0.020	26.08.97

8.2.3 Erasing PhotoCheck standards





At least 1 standard must still be stored to be able to perform the AQA check function (meter monitoring).

```
PhotoCheck-Standards
input
output
erases
return
```

- In the *PhotoCheck standards* submenu, call up the *erase* menu item.

```
erase PhotoCheck
445-2
445-3
445-4
return
```

Displays the stored PhotoCheck standards:

- Select with  or  
- Quit via *return*
- Erase with .


8. Analytical Quality Assurance (AQA)

8.2.4 Performing Photometer monitoring

Photometer monitoring (AQA1) includes a check of the

- Light barriers using the L1/L2 cells (contained within the scope of delivery of the Spectroquant® PhotoCheck)
- Absorbance measurement using PhotoCheck


standards.

- Press  to call up the *setup* menu
- Call up the *AQA check* submenu
- Call up the *meter* submenu.

The following display appears:

```
L-Check
-----
use L1

↵ cancel
```

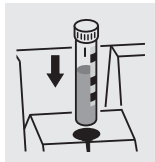

After approx.
1 s:

```
L-Check
-----
L1 ok
```

- Insert the L1 cell.




If the *error* message appears, clean the cell shaft with a damp, lint-free cloth and repeat the check.
If the message reoccurs, inform the service department.

After approx.
4 s:

```
L-Check
-----
use L2

↵ cancel
```

- Insert the L2 cell.

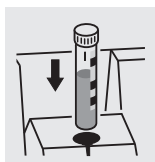


After approx.
1 s:

```
L-Check
-----
L2 ok
```

8. Analytical Quality Assurance (AQA)

After successful light barrier testing, the PhotoCheck standards (test solutions) are measured.

Example:

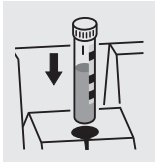

```
PhotoCheck      445-1
-----
use test
solution      445-1
↵ cancel
```

- Insert a cell with the test solution, 445-1. The photometer measures the absorbance of the test solution and compares the result with the value entered.

Absorbance test OK...

After approx. 3 s:

```
PhotoCheck      445-1
-----
                0.211 A
                ok
```

...or error message

```
PhotoCheck      445-1
-----
                ---- A
                Error
```

- Insert the next test solution
- Cancel:
To cancel the check means no release for the next "meter" AQA interval!

Error elimination:

1. Repeat the measurement (insert the cell again)
2. If necessary, perform a zero adjustment and repeat the check
3. Exchange the test solution (each packet contains two identical test solutions)
4. Use a new Spectroquant® Photo-Check packet
5. Quit and have the photometer checked in the factory

The absorbance test is terminated if an error message occurs and the meter is **not released**. On switching on, the warning message "AQA interval expired" appears until the AQA was successfully performed or the AQA mode was switched off.

Example: Report output

AQA check meter			AQA1		
26.08.97			10:23		
operator:					
AQA interval			12 weeks		
AQA check AQA1			ok		
L check			ok		
test sol.	meas. value	unit	theor. val.	tolerance	result
445-1	0.211	A	0.200	0.020	ok

8. Analytical Quality Assurance (AQA)

8.3 Total system monitoring with standard solutions (AQA2)

8.3.1 Entering standards



The standards compiled in the table "Spectroquant® CombiCheck and standard solutions" (see part "General information") are already stored method-specifically in the photometer. These values can be overwritten.

For **total system monitoring** (AQA2), only one standard per test can be stored at a time. The input of a standard is only complete with the input of the tolerances for finding it again, i.e. it is then first stored (no premature quitting).

- Press to call up the *setup* menu
- Call up the *meter setup* submenu.
- Call up the *AQA functions* submenu
- Input the password
- Call up the *AQA standards* submenu and the following display appears:

```
AQA-Standards
PhotoCheck
▶standard solution
return
```

- Call up the *standard solutions* submenu.

```
standard solution
▶input
output
erase
return
```

Select between

- *input*
Enter standards
- *output*
Print/display standards
- *erase*
Erase standards.

```
input standard

method: 2386
                14729
                P04-P
↓          0.5-25.0 mg/l
```

Displays the last selected method.

- Select the method with or
- Confirm with
- Input the standards.


8. Analytical Quality Assurance (AQA)

Example:

Method 14729 with a preset theoretical value of 15.0 mg/l and tolerance of 1.0 mg/l (CombiCheck 80).


Change to: theoretical value = 8 mg/l, tolerance = 0.7 mg/l (CombiCheck 20).

```
input standard
-----
method: 186
                14729
                PO4-P
↵      0.5-25.0 mg/l
```

- Confirm with .


```
standard      14729
-----
theor.val.: 15.0 mg/l
             (06.3-18.8 mg/l)
↵ confirm
```

- Enter the new theoretical value, e.g. 8.0 mg/l, via the numeric keypad.
Values in parentheses indicate the range in which the theoretical value should move.

- Confirm with .


```
standard      14729
-----
theor.val.: 08.0 mg/l
Tolerance: ±1.0 mg/l
↵ confirm
```

- Input the tolerance (0.7 mg/l) via the numeric keypad.

- Confirm with .

```
standard      14729
-----
theor.val.: 08.0 mg/l
Tolerance: ±0.7 mg/l
↵ confirm
```

Both standard and tolerance values have been overwritten.

- Confirm with .

8. Analytical Quality Assurance (AQA)

8.3.2 Output of standards

The current list of stored standards is output via the RS 232 interface (PC/printer) or via the display.

```
standard solution
input
output
erase
return
```

– Select the *output* submenu

– Confirm with .

```
download Standard
to display
to printer/PC
return
```

Select the output medium:

● to *display*

● to *printer/PC* (serial interface).

– Select with  or  

– Confirm with  to start the download.

Example: Report output

AQA check system		AQA2		
26.08.97		13:57		
system locked		on		
method	unit	theor. val.	tolerance	AQA date
14554	mg/l	2.00	0.20	24.08.97
14555	mg/l	5000	400	26.08.97

8. Analytical Quality Assurance (AQA)



8.3.3 Erasing standards

Erasing the method-specific standard solutions leads to the change of the measured value identification from AQA2 to AQA1 (with activated AQA mode).





```
AQA-Standards
-----
PhotoCheck
▶standard solution
return
```

- Call up the *standard solutions* submenu.

```
standard solution
-----
input
output
▶erase
return
```

- Select the menu item, *erase* with 
- Confirm with 

```
erase standard
-----
▶14560
14729
return
```


- Select the standard to be erased with  or  
- Erase with 

8.3.4 Monitoring of the total system using standard solutions

The AQA2 check can be performed after it has been activated (see section 8.1).
The following display appears:

```
AQA-Check
-----
insert cell
or
◀ start measurement
```

i For AQA2 with the setting, *n measurements*, we recommend to check and, if necessary, change the currently set number of *measurements* before each AQA check (8.1.3 CHANGING AQA INTERVALS).

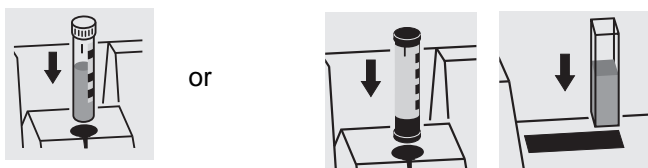

- Insert cell with prepared solution ready to be measured (e.g. using Spectroquant® CombiCheck). The photometer reads the barcode, identifies the method and performs the AQA2 check.

8. Analytical Quality Assurance (AQA)

After approx.
2 s:

AQA check OK ...

AQA-Check	14554
1.93 mg/l	
ok	

...or error message

AQA-Check	14554
3.45 mg/l	
Error	

- Repeat the check
If the error is repeated, perform troubleshooting of the error. See "Analytical Quality Assurance" in part "General information".



The *system* AQA2 check must be performed **separately** for each method monitored. The release is stored with the date and the specified interval. The AQA2 interval *system* set up for the respective method begins again.

Example: Report output (AQA mode: *n* weeks)

AQA check system			AQA2		
26.08.97			11:02		
operator:					
AQA interval			4 weeks		
method	meas. value	unit	theor. val.	tolerance	result
14554	1.95	mg/l	2.00	0.20	ok

Example: Report output (AQA mode: *n* measurements)


AQA check system			AQA2		
26.08.97			11:02		
operator:					
AQA interval			100 measurements		
method	meas. value	unit	theor. val.	tolerance	result
14554	1.95	mg/l	2.00	0.20	ok

9. Kinetics

The *kinetics* function enables the tracking of time-dependent changes in concentration or absorbance (for one selected wavelength) by repeating measurements over an adjustable time interval (setting possibility 00:05 to 60:00 min, smallest interval 5 seconds).

The kinetic measured values are stored in the memory in the same way as the measured values from single measurements (older measured values are overwritten). The maximum number of measuring cycles of a kinetic measurement is 1000 measurements (until the memory is full of kinetic measured values). The current measured value number appears in the header of the display. When the memory is full, the meter aborts the kinetic measurement. The starting values of the kinetic are not overwritten!

The output of kinetic measured values after the termination of the measurement is performed via the menu item, *download memory*.

- To switch on the photometer, open the cover.
- Press .
- In the *setup* menu, call up the *kinetics* submenu. The following display appears:

```
Kinetic
-----
Concentration
Absorbance
▶set interval
return
```


```
meas.interval
-----
set interval

  01 min 00 s

← confirm
```

```
Kinetic
-----
Concentration
Absorbance
▶set interval
return
```


- Call up the *set interval* menu item.

- Input the interval via the numeric keypad (default: 1 min).
Possible settings: 00:05 to 60:00
- Confirm with .

Call up the required measuring mode:

- Perform concentration measurement according to chapter 4
- Perform absorbance measurement according to chapter 5 .

10. Correction functions

- To switch on the photometer, open the cover.
- Press .
- In the *setup* menu, call up the *meter setup* submenu.
The following display appears:

```
Meter Setup
-----
return
AQA Functions
▶Correction Funct.
adjust zero
user def. methods
```

- Call up the *correction funct.* submenu.

The following display appears:

```
Correction Funct.
-----
▶Blank Value
Reference Absorbance
Turbidity Correct.
return
```

Select the correction function:

- blank value
- reference absorbance
- turbidity correct.

- Confirm with .

10.1 Blank value

The blank value (= reagent blank value) for each method is stored in the photometer. When the *blank value* function is active, the stored value is ignored and the measured value of a self-prepared reagent blank solution is used instead.

This procedure increases the measuring accuracy for some tests (for more information, see part "Analytical procedures").

A blank value is always stored for the method that was just called up. A maximum of 10 measured blank values can be stored, each of which is permanently assigned to a method.

A blank value remains stored until it is erased (menu item, *erase blank value*) or overwritten.


The *reset setup* function sets the *blank value* to *off*. The stored blank values, however, remain stored. The *reset total* function resets all settings and blank values at once.


If a measured blank value is stored and the *blank value* function is active for a method, this blank value is used for determining the measured value and the measured value is documented accordingly.

The *blank value* function is not active when delivered.

10. Correction functions

Measuring the concentration with a blank value


- Press  to call up the *concentration* measuring mode.

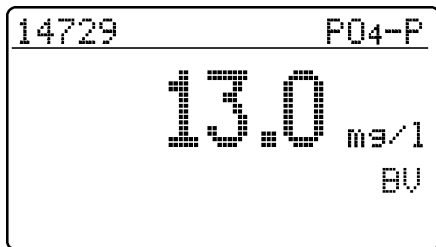





The value measured against the prepared blank solution is displayed.

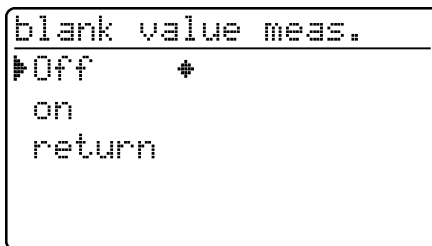10.1.1 Activating the blank value measurement

- In the *correction funct.* menu, call up the *blank value* submenu. The following display appears:




The *blank value meas.* function appears:

- Select the *on* menu item with  or  
- Confirm with .

i The stored blank values determined from blank solutions prepared by the user can be deactivated by switching off the blank value measurement. When doing so, the blank values remain stored in the memory and can be reactivated later.

Activating or deactivating the blank value function applies to all measurements using methods a blank value was stored for in the memory.

The *blank value* function is active and appears in the *setup* menu:

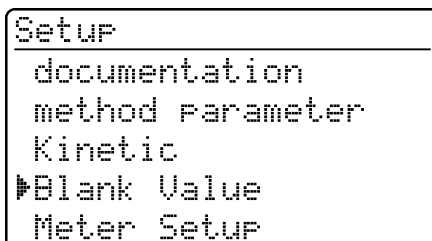

- To measure the blank value, call up the *blank value* submenu in the *setup* menu.

10. Correction functions

10.1.2 Measuring the blank value

```
Blank Value
└─ meas. blank value
   erase blank value
   recall blank values
   return
```

- Call up the *meas. blank value* menu item.



The menu items, *erase blank value* and *recall blank values* first appear after at least one blank value has been measured.

```
meas. blank value
      insert cell
      or
      ↵ start measurement
```

- Insert a cell with blank solution to start a measurement. The message, *measuring...*, appears on the display.

After approx.
2 s:

```
Blank Value      14729
      0.033 A
      ↵ return
```

Or, if all 10 storage locations for blank values are already occupied:

```
meas. blank value
      error: blank value
      memory full
      ↵ return
```



If the blank value of a method for which a blank value was already stored is measured again, this error message does not appear and the new measured value replaces the old measured value (new date as well).

- Selecting the menu item, *return*, returns to the menu item, *blank value*
- Before measuring and storing erase the old blank value.

```
Blank Value
└─ meas. blank value
   ↵ erase blank value
   recall blank values
   return
```

10.1.3 Erasing blank values

A measured blank value is erased via the menu item, *erase blank value*.


10. Correction functions

```
Blank Value
meas. blank value
erase blank value
recall blank values
return
```

```
erase blank value
all
single
return
```

```
erase blank value
12.01.2004 14:57
14758
0.100 A
```

```
erase blank value
12.01.2004 14:57
14758
erase
cancel
```

- Select the *erase blank value* menu item
- After confirming with  the *erase blank value* menu opens.











The *erase blank value* menu item first appears after a blank value has been measured.

Select between

- *all*
Erase all stored blank values
- *single*
Erase individual stored blank value







Each stored blank value is displayed with the date of the blank value measurement and the relevant method designation.





- Select the blank value with  or  
- Erase the displayed blank value with .
- Select the *erase* menu item with  or  
- Confirm with .

10.1.4 Recalling blank values

```
Blank Value
meas. blank value
erase blank value
recall blank values
return
```

```
recall blank values
12.01.2004 14:57
14758
0.100 A
return
```

- Select the *recall blank values* menu item with  or  
- Confirm with .

- Select the blank value with  or  
- Return with .

10. Correction functions

10.2 Reference absorbance

Each absorbance measurement is made against the basic absorbance stored in the meter. When activating the *reference absorbance*, this value is ignored and the value measured as the reference absorbance is used instead. The *reference absorbance* function is not active when delivered.

- In the *correction funct.* menu, call up the *reference absorbance* submenu.

The following display appears:

```
Reference Absorbance
▶Off      +
  on
  return
```

The *reference absorbance* function is active and appears in the *setup* menu:

```
Setup
-----
documentation
method parameter
Kinetic
▶ReferenceAbsorbance
Meter Setup
```

```
Reference Absorbance
▶measure ref. absorb.
  erase ref. absorb.
  return
```





```
Reference Absorbance

  insert cell
  or
  ↵ start measurement
```

The measured *reference absorbance value* remains stored until

- the photometer is switched off
- the method is changed
- it is manually erased using the *erase ref. absorb.* function.

The *reference absorbance* function appears:

- Select the *on* menu item with  or  
- Confirm with .

- To measure the reference absorbance, call up the *setup* submenu in the *reference absorbance* menu.

- Call up the *measure ref. absorb.* menu item.



The *erase ref. absorb.* menu item first appears after a reference absorbance has been measured.

- Insert the cell to start the measurement. The message, *measuring...*, ... appears on the display.


10. Correction functions

After approx.
2 s:

```

Ref. Absorbance ▶ 690nm
0.414 A
↵ return
    
```


The measured value for the wavelength displayed at the top right appears.

- If necessary, change the wavelength by pressing the  key.

The measured value for the selected wavelength appears.



The measured reference absorbance remains valid for all consecutive measurements with the same wavelength.

- Press  to call up the *absorbance* measuring mode.

```

meas. blank value
insert cell
or
↵ start measurement
    
```

- Insert the measuring cell.

After approx.
2 s:

```

Absorbance ▶ 690nm
0.345 A
Ref-A:0.414
    
```

The measured value corrected by the reference absorbance is displayed.


Erasing the reference value

A measured reference absorbance value is erased

- manually using the *erase ref. absorb.* function.
- by switching off the meter.

```

Reference Absorbance
measure ref. absorb.
▶ erase ref. absorb.
return
    
```

- Select the *erase ref. absorb.* menu item
- After confirming with , the reference absorbance value is erased.

10. Correction functions


10.3 Turbidity correction

Turbidity correction is used in sample solutions that contain finely distributed suspended particles. The suspended particles cause a light absorption. This leads to incorrect (too high) measured values. The function remains permanently switched on after it has been activated. Values that were measured using turbidity correction are given an identifier in the **display** and in the **documentation** (printout and storage).

- In the *correction funct.* menu, call up the *turbidity correct.* submenu.

The following display appears:

```
Turbidity Correct.
└─Off ─┬─
      on
      return
```

- Press  to call up the *concentration* measuring mode.

```
Concentration
-----
insert cell
or
start measurement
```

After approx.
2 s:

```
14542                               NO3-N
-----
          13.3 mg/l
TC
```

Warning of excessive turbidity:

If the turbidity absorbance of $0.100 A$ is exceeded, the meter displays the measured value together with a warning.





```
14542                               NO3-N
-----
          high turbidity!
TC          43 mg/l
```

The *turbidity correct.* function is not active when delivered.



This function is not necessary, or useful, in all methods. If the turbidity correction is active, the photometer automatically decides whether to perform the function or not depending on the method.

The *turbidity correct.* function appears:

- Select the *on* menu item with  or  
- Confirm with .

- Insert the measuring cell.

Display of the measured value with turbidity correction switched on: Identified by *TC*.

11. Zero adjustment

Zero adjustment is necessary

- after changing the lamp
- after the error message, *PhotoCheck* (AQA1) occurs
- on initial commissioning
- if the photometer was mechanically stressed, e.g. percussion, transport
- if the ambient temperature changed by more than 5 °C since the last zero adjustment
- at least every six months.

When performing the zero adjustment with a **round cell** observe the following points:

- Only use a clean, scratch-free round cell with distilled water. A prepared zero cell is provided with your photometer. In addition, a prepared zero cell is contained in the scope of delivery of the *PhotoCheck* (article 14693).
- If the round cell is visibly contaminated, or at least every 24 months, clean and refill it (minimum filling level 20 mm). Then check the cell for scratches.

When performing the zero adjustment with a **rectangular cell** observe the following points:

- With rectangular cells, the adjustment must be made using the same cell type (manufacturer) as for the measurement. This is important because glass of different manufacturers has different absorption behavior. When changing the cell type always repeat the zero adjustment with the new type.
- Clean the rectangular cell before performing the zero adjustment and fill it with distilled water (minimum filling level 20 mm).
- For measurement, always insert rectangular cells in the cell shaft with the same orientation as during the zero adjustment (e.g. cell labeling always on the left side).

11. Zero adjustment



Only perform the zero adjustment against distilled water in an optically perfect cell.

- Press
- In the *setup* menu, call up the *meter setup* submenu.
The following display appears:

```
Meter Setup
-----
reset
select language
set date/time
▶ adjust zero
user def. methods
```

- Call up the *zero adjustment* submenu with or .

```
adjust zero
-----
insert cell
or
↵ start measurement
```

- Insert a cell with distilled water.
The message, *measuring...*, appears on the display.

After approx.
2 s:

```
adjust zero
-----
10 mm ok
```

Successful zero adjustment for the 10 mm rectangular cell.




The zero adjustment must be performed separately for each cell type.

12. User-defined methods

User-defined methods are stored under code numbers. Numbers from 301 to 399 are allowed. These code numbers are used to fast-find user-defined methods when setting methods. A total of 50 user-defined methods can be stored. The 51st user-defined method causes the *method error* message to appear; if this occurs, erase an old method.

If you input a method designation that is already stored in the photometer, the characteristics are displayed and you can modify them. If the input is successful, the meter displays the accepted method.


i The photometer automatically restricts the measuring range (no message) if the entered method data allows absorbance values > 3.2 A.



- Press 
- In the *setup* menu, call up the *meter setup* submenu.
The following display appears:

```
Meter Setup
-----
return
AQA Functions
Correction Funct.
adjust zero
▶user def. methods
```

```
user def. methods
-----
▶input charact.
print charact.
erase
return
```


Measuring with user-defined methods:

```
select method
-----
method: 
TEST1
C6H5OH
0.1-22.3 mmol/l
```

- Call up the *user def. methods* submenu with  or .

Submenu *user def. methods* appears:




- *input charact.*
Entering characteristics for user-defined methods.
- *print charact.*
Printing the characteristics for user-defined methods; first appears after characteristics have been input.
- *erase*
Erasing single or all user-defined methods; first appears after characteristics have been input.

- Insert the measuring cell
- Input a method number via the numeric keypad
- Confirming with  starts the measurement.

12. User-defined methods


12.1 Entering characteristics via the keyboard

```
user def. methods
└─ input charact.
   print charact.
   erase
   return
```

- Select the *input charact.* submenu with  or  .



```
Input charact.
method number:
```

```
█_
```

- Input the method number (301 to 399) via the numeric keypad
- Confirm with . Confirming without character input results in a blank.

```
input characteristics
method-designation:
```



```
█_ _ _ _
```

- Input the method designation, e.g. TEST1:
 - Enter letters using , numerals via the numeric keypad
 - Confirm each time with .



Example:

```
input charact.
method designation:
█TEST1
```

```
input charact.
Wavelength:
 620 nm
 665 nm
└─ 690 nm
```



- Select the wavelength with 
- Confirm with .

```
input charact.
Citation:
C6H5OH
```

- Enter the citation form:
- Enter letters using , numerals via the numeric keypad
- Confirm each time with .

12. User-defined methods

```
input charact.  
-----  
unit:  
mmol/l
```



- Input the unit with 
- Confirm each time with 

```
input charact.  
-----  
Resolution:  
▶0.1  
  1  
 10
```



Selecting the resolution defines the numeric display of the beginning and end of the measuring range.

Possible selections:



- 0.001
- 0.01
- 0.1
- 1
- 10
- 100

- Select the resolution using 
- Confirm with 


```
input charact.  
-----  
zero point (E0):  
+0.009
```

- Select the sign (+/-) using 
- Enter the zero point:
 - Numbers from 0 to 9
 - Confirm each time with 


```
input charact.  
-----  
slope (E/C):  
+2.12
```

- Select the sign (+/-) using 
- Enter the slope:
 - Numbers from 0 to 9
 - Confirm each time with 

```
input charact.  
-----  
meas.range begin:  
  0.1  
mmol/l
```



- Enter the meas. range begin:
 - Numbers from 0 to 9
 - Confirm each time with 

```
input charact.  
-----  
meas.range end:  
 22.3  
mmol/l
```

- Enter the meas. range end:
 - Numbers from 0 to 9
 - Confirm each time with 

12. User-defined methods



```
input charact.
Reference Cell:
 50 mm
round
10 mm
```

- round
- 10 mm
- 20 mm
- 50 mm
- Select with 
- Confirm with 

```
method stored
Code.No.      301
designation    TEST1
Wavelength    690nm
unit          mmol/l
Citation      C6H5OH
```

Example:

Display and automatic printout of code no. and characteristics as a list.

- Scroll with 
- Return with 

12.2 Entering characteristics via the PC

Data format of the method data:

The data of user-defined methods is transferred in a string. The individual datablocks of the string are separated by spaces:



Further information is given in the chapter 15 RS 232 C INTERFACE.

Datablock	Character	Example
input user-def. methods	4	CEME
method number	3	301
method designation	5	TEST1
wavelength	5	690nm
unit	9	mmol/l
citation	12	C6H5OH
zero point	5	0.009
slope	5	2.12
meas. range begin	5	0.1
meas. range end	5	22.3
reference cell	2	(=round)
resolution	5	0.1

Example:

```
CEME 301 Test1 690nm mmol/l C6H5OH 0.009 2.12 0.1 22.3 14 0.1 <CR>
```

Error message "INVALID COMMAND"

The following inputs lead to this error message:

- slope = 0 or < -32000, > 32000
- zero point > 32000
- meas. range begin < 0 or > 32000
- meas. range end < 0 or > 32000
- meas. range end ≤ meas. range begin
- reference cell not equal to 10 mm, 14 mm, 20 mm or 50 mm
- designation of the wavelength does not agree with the photometer
- resolution not 0.1, 0.01, ...


12. User-defined methods

12.3 Printing characteristics

```
user def. methods
input charact.
print charact.
erase
return
```

– Call up the *print charact.* submenu with  or




– Start the print with .

```
print charact.
printing...
3 of 10 methods
printed

cancel
```

The characteristics of all user-defined methods are printed sequentially as a list.




– Cancel with .

Example: Report printout

< date >	< time >
user def. methods:	
code. no.	301
designation	TEST1
wavelength	690 nm
unit	mmol/l
citation	C6H5OH
resolution	0.1
zero point	+0.009
slope	+2.12
MRB	0.1
MRE	22.3
reference cell	round

12.4 Erasing user-defined methods


```
user def. methods
input charact.
print charact.
erase
return
```

– Call up the *erase* submenu with  or  .

12. User-defined methods





```
user def. methods
└─all
  single
  return
```

Select the required menu item:

- *all*
Erase all user-defined methods
 - *single*
Selectively erase individual methods
- Confirm with .



Selecting "all"

```
user def. methods
  all
└─erase
  cancel
```





- Call up the *erase* menu item with  or  
 - Confirm with .
- The message *methods are erased ...* appears

Selecting "single"


```
user def. methods
method: 301
                TEST1
                C6H5OH
└─ 1.0-24.2 mg/l
```

- Select a method, e.g. TEST1:
- Input the method number (301) via the numeric keypad
or
- by scrolling with 
- Confirm with .

```
user def. methods
  method
  301
└─erase
  cancel
```

- Call up the *erase* menu item with  or  
 - Confirm with .
- The message *method is erased ...* appears

13. Meter Setup

- To switch on the photometer, open the cover.
- Press 
- In the *setup* menu, call up the *meter setup* submenu.
The following display appears:

```
Meter Setup
-----
return
▶ AQA Functions
Correction Funct.
adjust zero
user def. methods
```

This chapter describes four functions of the *meter setup* menu:

- *select language*
- *set date/time*
- *Performing a meter reset*
- *system info*

13.1 Selecting the language

The following languages are stored in the photometer:

- Deutsch (German)
- English
- Français (French)
- Italiano (Italian)
- Português (Portuguese)
- Polski (Polish)
- Dansk (Danish)
- Svenska (Swedish)
- Español (Spanish)
- Nederlands (Dutch)
- Indonesia (Indonesian)
- Čeština (Czech)
- Magyar (Hungarian)
- Russkij (Russian)
- Türkçe (Turkish)
- Brasil (Brazilian)



This is the order in which the available languages appear in the *select language* menu.



The available languages are listed in the language of the respective country in the photometer.

When *Russkij* is selected as the language, the Cyrillic alphabet is used for the user guidance. Method designation and ID numbers are always displayed in Latin script. For output to the RS 232 C interface, Cyrillic characters are converted to Latin characters according to GOST.

```
meter setup
-----
adjust zero
user def. methods
set date/time
▶select language
system info
```

- Call up the *select language* menu item.

```
select language
-----
return
Deutsch
▶English
Français
Italiano
```

- Select a language, e.g. English
- Confirm with 
- Press the  key again:
Return to the *meter setup* submenu.
The displays appear in English.



13. Meter Setup

13.2 Setting the date/time

```
meter setup
-----
correction funct.
adjust zero
user def. methods
▶set date/time
select language
```

- Call up the *set date/time* menu item.

```
Date/Time
-----
Date      01.01.98
          (dd.mm.yy)
Time      16:45
          (hh:mm)
⏏ confirm
```

- Input the date via the numeric keypad
- Confirm with 
- Input the time via the numeric keypad
- Confirm with 

13.3 Reset

It is possible to reset the photometer to its factory settings (delivery state) in single steps. The *reset total* function resets all settings and blank values at once.



All AQA functions are retained when *meter setup* is used. See section 8.1.6 for AQA reset.

```
Meter Setup
-----
set date/time
select language
system info
▶reset
return
```

- Call up the *reset* menu item.

```
reset
-----
▶total
meas.storage
Setup
return
```

Select between

- *total*
Erase the measured value storage and reset the settings to the delivery state
- *meas. storage*
Erase the measured value storage
- *setup*
Reset all settings to the delivery state.

Example: Performing a total reset

```
reset
-----
reset total
▶reset
cancel
```

- Select the *reset* menu item

- Confirm with 

A meter reset is performed (measured value memory and setup).

13. Meter Setup

13.4 System info

```
Meter Setup
-----
user def. methods
set date/time
select language
▶system info
reset
```

- Call up the *system info* menu item.

```
Meter Setup
-----

Software: 2.01
methods: 19.00

↵ return
```

Sample display

14. Updating method data

You will always find the latest method data for your photometer on the Internet. A method update contains all new test sets and methods respectively. Additionally, minor modifications of already existing methods are transferred with it. With a method update, you receive all new methods and, at the same time, can easily and conveniently update all method data.

The software provided for downloading contains the program file and method data. It can be downloaded from our homepage with a mouse click.

The files are packed in a self-decompressing archive file (*.exe) or in a zip file (*.zip) and can be decompressed after the download.

Carry out the update as follows:

To download and update the photometer method data via the built-in RS232 interface, you need the following:

- PC (Win 95 or higher) with Internet connection
- PC cable (available as an accessory)
- An *.exe or *.zip file from the Internet; contains the "UpdateMethodData.exe" program file and 6 method data files (pls6md.xxx, pls12md.xxx, plspekmd.xxx, nova30md.xxx, nova60md.xxx, nova400md.xxx; xxx = version).

- Switch on the photometer (open the cover).
- Switch on the PC.
- Download from the Internet the software including the method data (*.exe or *.zip) and copy it into a separate directory or on a floppy disk.
- Decompress the *.exe file with a double-click or decompress the *.zip file with Winzip.
- Connect the serial interfaces of the PC and photometer with the cable.
- Start the "UpdateMethodData.exe" program file by double-clicking. The "Update Method Data" window appears. In the upper half of the window there is the name of your photometer (among other things), behind it there is the method version in brackets (e.g. 8.00).



All method data are reloaded into the photometer with the update. The old method data are overwritten by this.

- Click on the "Search meter" button. The program automatically recognizes the connected photometer. Another "Update Method Data" window appears.
- Click on the "Start" button to start the method download. The process takes approx. 3 minutes. You can terminate it at any time by clicking on the "Cancel" button. In this case, however, the download has to be carried out once again completely so that the photometer can save the method data and is operative.

14. Updating method data

During the download, the following display appears on the photometer screen:

```
remote
```

- After the download, confirm the "Data successfully downloaded" message. The download is finished. The photometer returns to the *concentration* measuring mode.



You can check whether the new method data are stored in the photometer. To do so, proceed as follows:

```
Meter Setup
-----
user def. methods
set date/time
select language
▶system info
reset
```

```
Meter Setup
-----
Software: 2.01
methods: 19.00
↵ return
```

- In the setup menu, call up the meter setup submenu.
- Call up the *system info* menu item.

Sample display (the software version is irrelevant here). The method version (here: 19.00) has to agree with the method version for your photometer in the "Update Method Data" window during the download.

Error messages

Message	Meaning	Remedy
No meter found	Connection PC - photometer out of order or not available	<ul style="list-style-type: none">- Tightly connect the cable to the serial interfaces of the PC and photometer.- Use the correct cable
	Photometer not recognized	<ul style="list-style-type: none">- Select the photometer manually

15. RS 232 C interface

Via the interface, data can be

- output to a printer and
- exchanged with a personal computer (PC)

For this, the following items are available as accessories:

- Printer cable
- Printer
- Interface cable
- Communication software.

15.1 Principle course of the remote control

String to meter	Reply from meter	Operating mode
S <CR>	> <CR>	Remote (remote control)
Command xx (see 15.2 command list)	Reply string command xx <CR>	Remote (remote control)
.		
.		
.		
CLOC <CR>		Concentration measurement



The keyboard of the photometer is locked in the *remote* operating mode.

15.2 Command list

Command	Function
S	Begin communication
CLOC	Switchover to normal operation (concentration measurement)
CDAT [anz]	Reads out stored measured values; [anz] = number of the measured values to be output
CMES [MMM]	Measurement and transmission of the concentration value with date/time; [MMM] = method number (e.g. 086 for method 14729)
CEXT [LLL]	Measurement and transmission of the absorbance value for the wavelength; [LLL] = wavelength
CBLA [MMM]	Measurement and transmission of the sample blank value; [MMM] = method number
CCLB [MMM]	Erase measured sample blank values; [MMM] = method number
CEME	Input user-defined methods (see chapter chapter 12 USER-DEFINED METHODS)
REME [MMM]	Output user-defined methods; [MMM] = method number



The error message, *Invalid command*, appears if commands are unknown or cannot be carried out (e. g. if optional parameters do not agree with the cell coding). Optional parameters [MMM] and [LLL] need only be input for uncoded cells.

15. RS 232 C interface

15.3 Output format of measured values

Character	Meaning
3	consecutive number (not required for interface commands CMES, CEXT and CBLA)
5	method designation
6	I. D. number
17	date and time
4	special characters
9	meas. value
10	unit
12	citation
4	AQA ID (AQA2/AQA1)
4	dilution factor

Notes:

Data fields are separated by spaces.
Character set: IBM, code page 437


Meaning of the special characters:

- ! = Measuring with blank value (concentration) or reference absorbance (absorbance)
- t/T = Measurement with turbidity correction/with high turbidity
- * = Measured value outside the measuring range
- Q = AQA measurement

15.4 Data transmission

Baud rate	4800
Data bits:	8
Stop bits:	1
Parity:	none
Handshake:	Hardware
Max. cable length	15 m


15.5 Pin assignment

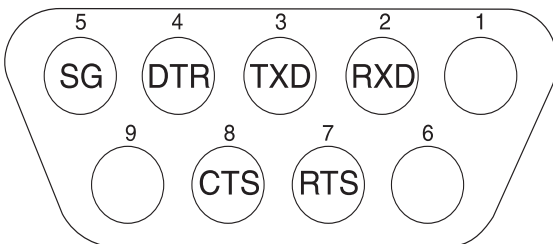

Photometer 9-pin socket	Computer 9-pin socket	25 pin plug	Printer with RS 232 C interface
1	4	20	-
2	3	2	TXD
3	2	3	RXD
4	1 and 6	6	-
5	5	7	SG
6	4	20	-
7	8	5	-
8	7	4	DTR (if not available: short-circuit CTS and RTS)
9		-	-

16. Maintenance

16.1 Changing the lamp

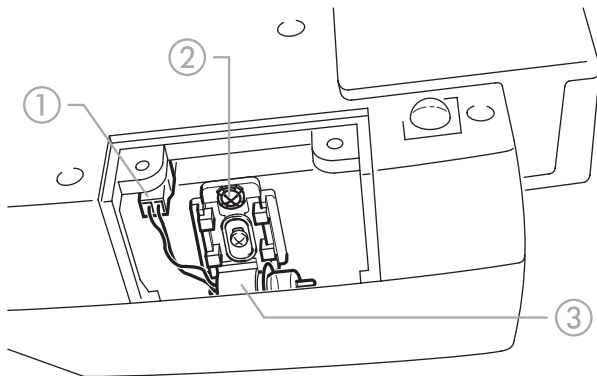

- Switch off the photometer and disconnect it from the power line
- Carefully turn up the photometer and park it safely
- Screw off the lamp cover on the underside of the photometer




Let the lamp of the photometer cool down.


- Pull out the plug ①
- Unscrew the screw ②
- Remove the lamp with its holder ③ by pulling it gently upwards



Do not touch the new light bulb of the photometer.

- Insert a new preset lamp and screw it tight using the screw ②

- Connect the plug ① of the new lamp
- Screw the lamp cover on again
- Set up the photometer again and connect it to the power line
- Press and hold 
- Switch on the meter (open the cover) and after the following display appears, release :

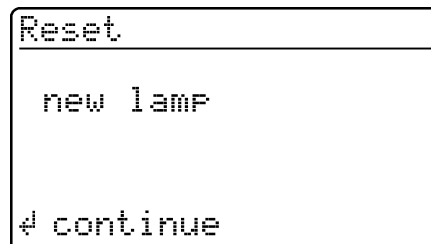

- Press 

16.2 Actions to take if a cell is broken



Do not rotate the photometer to pour out the liquid!

The photometer has a draining mechanism under the cell shaft that, when operated correctly, prevents any liquid coming into contact with electronic components.

- Switch off the photometer (close the cover) and disconnect it from the line power
- Let the liquid drain off
- Carefully remove any pieces of glass, e.g. using tweezers
- Carefully clean the cell shaft with a damp, lint-free cloth

- Let the cell shaft dry
After it is dry, check the photometer:
- Perform a photometer monitoring (see section 8.2).

17. Technical Data

Optical measuring principle	Filter photometer with reference beam absorption measurement; simultaneous recording of all wavelengths
Light source	Tungsten halogen lamp, preset
Receiver	12 x photo diode array
Optical filters	340 nm, 410 nm, 445 nm, 500 nm, 525 nm, 550 nm, 565 nm, 605 nm, 620 nm, 665 nm, 690 nm, 820 nm, Accuracy: ± 2 nm; Half width: 340 nm = 30 nm ± 2 nm; all others = 10 nm ± 2 nm
Photometric reproducibility	0.001 A at 1.000 A
Photometric resolution	0.001 A
Warm-up time	none
Measuring time	approx. 2 s
Types of measurement	Concentration (method dependent, selectable display form), absorbance, transmission
Measuring range absorbance	-0.300 A to 3.200 A
Measuring range transmission	0.1 % to 1000 %
Balancing	Permanently stored
Drift correction	Automatic on each Self-Check
Retrofitting of new methods	via the Internet
User-defined methods	max. 10
Kinetics	Automatic measurement repetition with selectable interval
Bar code recognition	automatic selection of the method; automatic recognition of the reagents lot
Cell recognition	automatic

Self-Check	<p><i>Test:</i> Memory, optics, electronic measured value recording, barcode recognition, cell recognition</p> <p><i>Automatic calibration:</i> Optics, electronic measured value recording, barcode recognition, rectangular cell recognition</p>
Time/Date	Real-time clock in the photometer
Dimensions	H: 140 mm, D: 270 mm, W: 260 mm
Weight	approx. 2.3 kg (battery version: 2.8 kg)
Meter safety	EN 61010, IEC 1010
Safety class	EN 61010-1/class 3
Power pack	
<ul style="list-style-type: none"> ● Type 	<p>Friwo FW6798/11.8363 * Friwo Part-No. 1810502 Input: 230 V~ $\pm 10\%$/50 Hz/25 VA Output: 12 V~/1540 mA</p> <p>Friwo FW6798/11.8365 * Friwo Part-No. 1769227 Input: 120 V~ $\pm 10\%$/60 Hz/24 VA Output: 12 V~/1540 mA</p> <p>* compulsory for meters with UL/cUL test certificates</p> <p>FRIWO FW 7555O/15 Friwo Part. No. 1822367 Input: 100 ... 240 V ~ / 50 ... 60 Hz / 400 mA Output: 15 V DC / 1 A</p>
<ul style="list-style-type: none"> ● Meter safety 	EN 60950
Battery operation (optional)	<p>Built-in battery: NiCad rechargeable battery 7.2 V/2200 mAh, operating time with new, fully charged battery: typical 40 hours with 10 measurements per hour, trickle charging in line operation, approx. 5 h charging time for a discharged battery, total discharge protection</p>
Power consumption in line operation	max. 1300 mA

17. Technical Data

EMC	EU directive 89/336/EEC EN 61326-1 EN 61000-3-2 A14 EN 61000-3-3 FCC class A
Climatic class	2, VDI/VDE 3540
Ambient temperature	Storage: -25 °C to +65 °C Operation: +5 °C to +40 °C
Allowable relative humidity	Annual mean: 75 % 30 days/year: 95 % other days: 85 %
Test certificate	CE
Operating elements	On/off switch actuated by opening/closing the lid of the cell shaft cover Silicon keyboard with 4 function keys and number field with 12 keys Round cell shaft – for round cells (flat cell floor, external/internal diameter 16 mm / 13.8 mm) Rectangular cell shaft – for 10 mm, 20 mm and 50 mm rectangular cells with a maximum width of 12.60 mm Recess for MemoChip
Display	Graphical display 128 x 64 pixels
Connections	
● Digital interface	RS 232 C 9-pin socket to connect to PC or printer
● Power supply	2-pin socket to connect the plug-in power supply unit
Data storage	Cyclical memory to record 1000 measured values

FCC Class A Equipment Statement

Note: This equipment has been tested and found to comply with the limits for a Class A digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference when the equipment is operated in a commercial environment. This equipment generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instruction manual, may cause harmful interference to radio communications. Operation of this equipment in a residential area is likely to cause harmful interference in which case the user will be required to correct the interference at his own expense.

17. Technical Data

Software settings when delivered

ID number input:	Off
Measured value number:	1
blank value is:	Off
Reference absorbance:	Off
turbidity correct.:	Off
language:	English
Kinetic interval:	60s
Date of the last valid AQA1 check:	invalid (not yet measured)
AQA1 interval:	12 weeks
AQA2 interval:	4 weeks
AQA password:	0000
AQA mode:	Off
Lock measurement if AQA2 expired:	Off
Checks to be measured with AQA1:	none
AQA2 values:	none

Settings after reset - total

Measured value storage and setup reset

Settings after reset - meas. storage

Meas. value number:	1
Measured values:	none


Settings after reset - setup

ID number input:	Off
Measured value number:	1
blank value:	Off
reference absorbance:	Off
turbidity correct.:	Off
Language:	unchanged
Kinetic interval:	60s
Citation forms of the methods	respective reference citation
Units of the methods:	respective reference unit

Settings after reset - AQA

Date of the last valid AQA1 check:	invalid (not yet measured)
AQA1 interval:	12 weeks
AQA2 interval:	4 weeks
AQA password:	0000
AQA mode:	Off
Lock measurement if AQA2 expired:	Off
Checks to be measured with AQA1:	none (Input theoretical values and tolerances are not erased and are offered again with the next input).
AQA2 values:	none (theoretical values and tolerances of all methods are set to default values according to the "Spectroquant® CombiCheck and standard solutions" table in the part "General information".)

18. What to do if...

The display remains blank when switched on	Connect the photometer to the line power via the power pack. In the case of battery operation: Battery empty, charging required (approx. 5h); line operation is possible without restrictions during charging time.
 appears	Battery nearly empty. Charging required (see chapter 3 COMMISSIONING).
Date/time is lost when switched off	The backup battery of the real time clock is empty and has to be replaced. Send the photometer to the service department for this.
MemoChip is not recognized	The MemoChip is not recognized by the photometer though it is in the recess during switching on. Operate the photometer with line power (see chapter 3 COMMISSIONING). Repeat the procedure.
Password forgotten	Inform the service department.
Photometer does not react	The connected printer is off line. Switch on the printer or pull out the interface cable.
Error messages:	
<i>remove cell</i>	The message remove cell appears on the display although no cell is inserted. Clean the cell shaft with a damp, lint-free cloth. If the error message still appears, return the photometer to the service department.
<i>lamp defective</i>	Replace the lamp (see chapter chapter 16 MAINTENANCE).
<i>no zero adjustment</i>	No zero adjustment is stored in the meter for the cell. Perform zero adjustment (see chapter chapter 11 ZERO ADJUSTMENT).
<i>cell error</i>	The rectangular cell is inserted incorrectly, or two cells are inserted. Insert the cell correctly.
<i>cell invalid</i>	A cell type that is not allowed for the selected method was selected, e. g. a round cell for the reagent test.
<i>method invalid</i>	No data is stored in the photometer for the selected method. Update method data (see chapter chapter 14 UPDATING METHOD DATA).
<i>wrong method</i>	During a difference measurement, the method was changed between the first and second measurement. During a difference measurement, the method must remain identical.
<i>E_0</i>	Hardware error: Send the photometer to the service department.
<i>E_1, E_2 or E_3</i>	Replace the lamp (see chapter chapter 16 MAINTENANCE). If the error message remains, send the meter to the service department.



Wissenschaftlich-Technische Werkstätten GmbH

Dr.-Karl-Slevogt-Straße 1
D-82362 Weilheim

Germany

Tel: +49 (0) 881 183-0
+49 (0) 881 183-100
Fax: +49 (0) 881 183-420
E-Mail: Info@WTW.com
Internet: <http://www.WTW.com>